

ACR Guidance Document for Safe MR Practices: 2007

Emanuel Kanal¹
 A. James Barkovich²
 Charlotte Bell³
 James P. Borgstede⁴
 William G. Bradley, Jr.⁵
 Jerry W. Froelich⁶
 Tobias Gilk⁷
 J. Rod Gimbel⁸
 John Gosbee⁹
 Ellisa Kuhni-Kaminski¹
 James W. Lester, Jr.¹⁰
 John Nyenhuis¹¹
 Yoav Parag¹
 Daniel J. Schaefer¹²
 Elizabeth A. Sebek-Scoumis¹
 Jeffrey Weinreb¹³
 Loren A. Zaremba¹⁴
 Pamela Wilcox¹⁵
 Leonard Lucey¹⁵
 Nancy Sass¹⁵
 for the ACR Blue Ribbon Panel on MR Safety

Keywords: MR contrast agents, MRI, safety

DOI:10.2214/AJR.06.1616

Received December 8, 2006; accepted without revision December 18, 2006.

E. Kanal is a consultant for, is a member of the speakers bureau of, and provides research support for Bracco Diagnostics and GE Healthcare; is a member of the speakers bureau of and provides research support for Siemens Medical Solutions; and provides research support for Berlex and Medtronic.
 T. Gilk is a consultant for Mednovus, Inc.
 J. R. Gimbel provides research support for St. Jude Medical, Medtronic, and Biotronik.
 J. Nyenhuis is a consultant for and provides research support to Medtronic.
 J. Weinreb is a consultant and member of the speakers bureau for GE Healthcare.

AJR 2007; 188:1447–1474

0361–803X/07/1886–1447

© American Roentgen Ray Society

There are potential risks in the MR environment, not only for the patient [1, 2] but also for the accompanying family members, attending health care professionals, and others who find themselves only occasionally or rarely in the magnetic fields of MR scanners, such as security or housekeeping personnel, firefighters, police, etc. [3–6]. There have been reports in the medical literature and print media detailing magnetic resonance imaging (MRI) adverse incidents involving patients, equipment, and personnel that spotlighted the need for a safety review by an expert panel. To this end, the American College of Radiology (ACR) originally formed the Blue Ribbon Panel on MR Safety. First constituted in 2001, the panel was charged with reviewing existing MR safe practices and guidelines [5–9] and issuing new ones as appropriate for MR examinations. Published initially in 2002 [3], the ACR MR Safe Practice Guidelines established de facto industry standards for safe and responsible practices in clinical and research MR environments. These were subsequently reviewed and updated in May 2004 [4]. After

reviewing substantial feedback from the field and installed bases, as well as changes that had transpired throughout the MR industry since the publication of the 2004 version of this document, the panel extensively reviewed, modified, and updated the entire document in 2006–2007.

The present panel consists of the following members: A. James Barkovich, MD; Charlotte Bell, MD (American Society of Anesthesiologists); James P. Borgstede, MD, FACC; William G. Bradley, MD, PhD, FACC; Jerry W. Froelich, MD; Tobias Gilk, architect; J. Rod Gimbel, MD, FACC, cardiologist; John Gosbee, MD, MS; Ellisa Kuhni-Kaminski, RT (R)(MR); Emanuel Kanal, MD, FACC, FISMRM (chair); James W. Lester, MD; John Nyenhuis, PhD; Yoav Parag, MD; Daniel Joe Schaefer, PhD, engineer; Elizabeth A. Sebek-Scoumis, RN, BSN, CRN; Jeffrey Weinreb, MD; Loren A. Zaremba, PhD, FDA; Pamela Wilcox, RN, MBA (ACR staff); Leonard Lucey, JD, LLM (ACR staff); and Nancy Sass, RT (R)(MR)(CT) (ACR staff). The following represents the most recently modified and updated version of the combined prior two re-

¹Department of Radiology, University of Pittsburgh Medical Center, Pittsburgh, PA.

²Neuroradiology Section, University of California San Francisco, San Francisco, CA.

³American Society of Anesthesiologists and Department of Anesthesiology, New York University School of Medicine, New York, NY.

⁴Colorado Springs Radiologists, Colorado Springs, CO.

⁵Professor and Chairman, Department of Radiology, University of California San Diego, San Diego, CA.

⁶Department of Radiology, University of Minnesota, Minneapolis, MN.

⁷MRI-Planning, Kansas City, MO.

⁸East Tennessee Heart Consultants, Lenoir City, TN.

⁹University of Michigan Health System and Red Forest Consulting, Ann Arbor, MI.

¹⁰Chapel Hill, NC.

¹¹Department of Electrical and Computer Engineering, Purdue University, West Lafayette, IN.

¹²MR Systems Engineering, GE Healthcare, Milwaukee, WI.

¹³Yale University School of Medicine, New Haven, CT.

¹⁴U.S. Food and Drug Administration, Rockville, MD.

¹⁵American College of Radiology, 1891 Preston White Dr., Reston, VA 20191. Address correspondence to N. Sass.

ports [3, 4] issued by the American College of Radiology Blue Ribbon Panel on MR Safety, chaired by Emanuel Kanal, MD, FACR. It is important to note that nothing that appears herein is the result of a “majority vote” of the members of this panel. As with each prior publication of these ACR MR Safe Practice Guidelines, the entire document, from introduction to the markedly expanded appendices, represents the unanimous consensus of each and every member of this Safety Committee and the various areas of expertise that they represent. This includes representation from fields and backgrounds as diverse as MR physicists, research/academic radiologists, private practice radiologists, MR safety experts, patient safety experts/researchers, MR technologists, MR nursing, National Electrical Manufacturers Association, the U.S. Food and Drug Administration (FDA), the American Society of Anesthesiologists, legal counsel, and others. Lay personnel, physicians, PhDs, department chairs and house-staff/residents, government employees and private practitioners, doctors, nurses, technologists, radiologists, anesthesiologists, cardiologists, attorneys—these are all represented on this Committee. It was felt that achieving unanimity for these guidelines was critical in order to demonstrate to all that these guidelines are not only appropriate from a scientific point of view, but are reasonably applicable in the real world in which we all must live, with all its patient care, financial, and throughput pressures and considerations.

The following MR safe practice guidelines document is intended to be used as a template for MR facilities to follow in the development of an MR safety program. These guidelines were developed to help guide MR practitioners regarding these issues and to provide a basis for them to develop and implement their own MR policies and practices. It is intended that these MR safe practice guidelines (and the policies and procedures to which they give rise) be reviewed and updated on a regular basis as the field of MR safety continues to evolve.

The principles behind these MR safe practice guidelines are specifically intended to apply not only to diagnostic settings but also to patient, research subject, and health care personnel safety for all MRI settings, including those designed for clinical diagnostic imaging, research, interventional, and intraoperative MR applications.

With the increasing advent and use of 3.0-Tesla and higher strength magnets, users need to recognize that one should never assume MR compatibility or safety information about a device if it is not clearly documented in writing. Decisions based on published MR safety and compatibility claims should recognize that all such claims apply only to specifically tested conditions, such as static magnetic field strengths, static gradient magnetic field strengths and spatial distributions, and the strengths and rates of change of gradient and radiofrequency (RF) magnetic fields.

Finally, there are many issues that impact MR safety that should be considered during site planning for a given MR installation. These have historically not been dealt with in the prior versions of the ACR MR Safe Practice Guidelines. For the first time, we include in this article, as separate appendices, sections that address such issues as well, including cryogen emergency vent locations and pathways, 5-gauss lines, siting considerations, patient access pathways, etc. Yet despite their appearance herein, these issues, and many others, should be reviewed with those experienced in MR site planning and familiar with the patient safety and patient flow considerations prior to committing to construction of a specific site design. In this regard, enlisting the assistance of an architectural firm experienced in this area, and doing so early in the design stages of the planning process, may prove most valuable.

It remains the intent of the ACR that these MR Safe Practice Guidelines will prove helpful as the field of MRI continues to evolve and mature, providing MR services that are among the most powerful, yet safest, of all diagnostic procedures to be developed in the history of modern medicine.

ACR Guidance Document for Safe MR Practices: 2007

A. Establish, Implement, and Maintain Current MR Safety Policies and Procedures

1. All clinical and research MR sites, irrespective of magnet format or field strength, including installations for diagnostic, research, interventional, and/or surgical applications, should maintain MR safety policies.
2. These policies and procedures should also be reviewed concurrently with the introduction of any significant changes in safety parameters of the MR environment of the site (e.g., adding faster or stronger gradient capabilities or higher RF duty cycle studies) and updated as needed. In this review process, national and international standards and recommendations should be taken into consideration prior to establishing local guidelines, policies, and procedures.
3. Each site will name an MR medical director whose responsibilities will include ensuring that MR safe practice guidelines are established and maintained as current and appropriate for the site. It is the responsibility of the site's administration to ensure that the policies and procedures that result from these MR safe practice guidelines are implemented and adhered to at all times by all of the site's personnel.
4. Procedures should be in place to ensure that any and all adverse events, MR safety incidents, or “near incidents” that occur in the MR site are reported to the medical director in a timely fashion (e.g., within 24 hours or 1 business day of their occurrence) and used in continuous quality improvement efforts. It should be stressed that the Food and Drug Administration states that it is incumbent upon the sites to also report adverse events and incidents to them via their MedWatch program. The ACR supports this requirement and feels that it is in the ultimate best interest of all MR practitioners to create and maintain this consolidated database of such events to help us all learn about them and how to better avoid them in the future [10, 11].

B. Static Magnetic Field Issues: Site Access Restriction

I. Zoning

The MR site is conceptually divided into four Zones (see Figure 1 and Appendix 1):

- a. Zone I: This region includes all areas that are freely accessible to the general public. This area is typically outside the MR environment itself and is the area through which patients, health care personnel, and other employees of the MR site access the MR environment.
- b. Zone II: This area is the interface between the publicly accessible, uncontrolled Zone I and the strictly controlled Zones III and IV. Typically, patients are greeted in Zone II and are not free to move throughout Zone II at will, but are rather under the supervision of MR personnel (see section B.2.b, below). It is in Zone II that the answers to MR screening questions, patient histories, medical insurance questions, etc. are typically obtained.

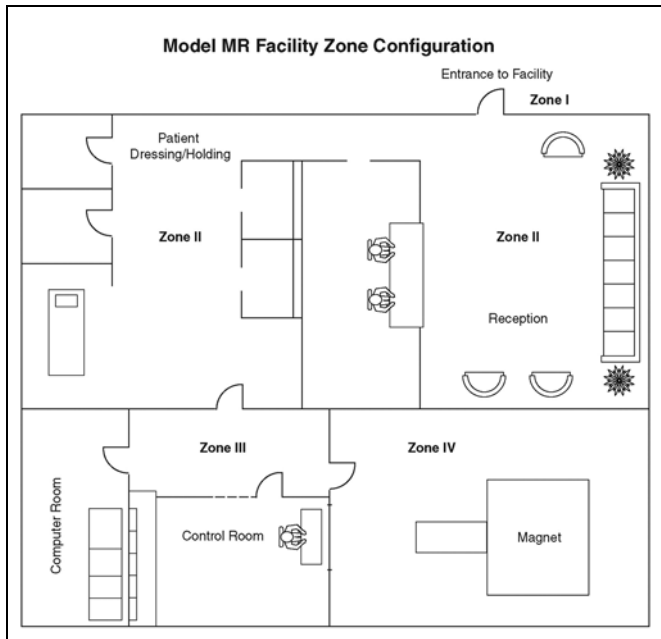


Fig. 1—Idealized sample floor plan illustrates site access restriction considerations. Other MR potential safety issues, such as magnet site planning related to fringe magnetic field considerations, are not meant to be included herein. See Appendix 1 for personnel and zone definitions. Note—In any zone of the facility, there should be compliance with Health Insurance Portability and Accountability Act (HIPAA) regulations in regard to privacy of patient information. However, in Zone III, there should be a privacy barrier so that unauthorized persons cannot view control panels.

- c. **Zone III:** This area is the region in which free access by unscreened non-MR personnel or ferromagnetic objects or equipment can result in serious injury or death as a result of interactions between the individuals or equipment and the MR scanner's particular environment. These interactions include, but are not limited to, those involving the MR scanner's static and time-varying magnetic fields. All access to Zone III is to be strictly restricted, with access to regions within it (including Zone IV, see below) controlled by, and entirely under the supervision of, MR personnel (see section B.2.b, below). Specifically identified MR personnel (typically, but not necessarily only, the MR technologists) are to be charged with ensuring that this MR safe practice guideline is strictly adhered to for the safety of the patients and other non-MR personnel, the health care personnel, and the equipment itself. This function of the MR personnel is directly under the authority and responsibility of the MR medical director or the level 2 MR personnel—designated (see section B.2.b, below) physician of the day for the MR site.

Zone III regions should be physically restricted from general public access by, for example, key locks, passkey locking systems, or any other reliable, physically restricting method that can differentiate between MR personnel and non-MR personnel. The use of combination locks is discouraged as combinations often become more widely distributed than initially intended, resulting in site restriction violations being more likely with these devices. Only MR personnel shall be provided free access, such as the access keys or passkeys, to Zone III.

There should be *no* exceptions to this guideline. Specifically, this includes hospital or site administration, physician, se-

curity, and other non-MR personnel (see section B.2.c, below). Non-MR personnel are not to be provided with independent Zone III access until such time as they undergo the proper education and training to become MR personnel themselves. Zone III, or at the very least the area within it wherein the static magnetic field's strength exceeds 5 gauss, should be demarcated and clearly marked as being potentially hazardous.

Because magnetic fields are three-dimensional volumes, Zone III controlled access areas may project through floors and ceilings of MRI suites, imposing magnetic field hazards on persons on floors other than that of the MR scanner. Zones of magnetic field hazard should be clearly delineated, even in typically nonoccupied areas such as rooftops or storage rooms, and access to these Zone III areas should be similarly restricted from non-MR personnel as they would be inside any other Zone III region associated with the MRI suite. For this reason, magnetic field strength plots for all MRI systems should be analyzed in vertical section as well as in horizontal plan, identifying areas above or below, in addition to areas on the same level, where persons may be at risk of interactions with the magnetic field.

- d. **Zone IV:** This area is synonymous with the MR scanner magnet room itself, that is, the physical confines of the room within which the MR scanner is located. Zone IV, by definition, will always be located within Zone III, as it is the MR magnet and its associated magnetic field that generates the existence of Zone III. Zone IV should also be demarcated and clearly marked as being potentially hazardous due to the presence of very strong magnetic fields. As part of the Zone IV site restriction, all MR installations should provide for direct visual observation by level 2 personnel to access pathways into Zone IV. By means of illustration only, the MR technologists would be able to directly observe and control, via line of sight or via video monitors, the entrances or access corridors to Zone IV from their normal positions when stationed at their desks in the scan control room.

Zone IV should be clearly marked with a red light and lighted sign stating, "The Magnet is On." Except for resistive systems, this light and sign should be illuminated at all times and should be provided with a backup energy source to continue to remain illuminated for at least 24 hours in the event of a loss of power to the site.

In case of cardiac or respiratory arrest or other medical emergency within Zone IV for which emergent medical intervention or resuscitation is required, appropriately trained and certified MR personnel should immediately initiate basic life support or CPR as required by the situation *while* the patient is being emergently removed from Zone IV to a predetermined, magnetically safe location. *All* priorities should be focused on stabilizing (e.g., basic life support with cardiac compressions and manual ventilation) and then evacuating the patient as rapidly and safely as possible from the magnetic environment that might restrict safe resuscitative efforts.

Further, for logistical safety reasons, the patient should always be moved from Zone IV to the prospectively identified location where full resuscitative efforts are to continue. (See Appendix 2.)

Quenching the magnet (for superconducting systems only) is not routinely advised for cardiac or respiratory arrest or other medical emergency, since quenching the magnet and having the magnetic field dissipate could easily take more than a minute. Further-

more, as quenching a magnet can theoretically be hazardous, ideally one should evacuate the magnet room, when possible, for an intentional quench. One should rather use that time wisely to initiate life support measures while removing the patient from Zone IV to a location where the strength of the magnetic field is insufficient to be a medical concern. Zones III and IV site access restriction *must* be maintained during resuscitation and other emergent situations for the protection of all involved.

2. MR personnel and non-MR personnel

- a. All individuals working within at least Zone III of the MR environment should be documented as having successfully completed at least one of the MR safety live lectures or prerecorded presentations approved by the MR medical director. Attendance should be repeated at least annually, and appropriate documentation should be provided to confirm these ongoing educational efforts. These individuals shall be referred to henceforth as MR personnel.
- b. There are two levels of MR personnel:
 1. Level 1 MR personnel: Those who have passed minimal safety educational efforts to ensure their own safety as they work within Zone III will be referred to henceforth as level 1 MR personnel.
 2. Level 2 MR personnel: Those who have been more extensively trained and educated in the broader aspects of MR safety issues, including, for example, issues related to the potential for thermal loading or burns and direct neuromuscular excitation from rapidly changing gradients, will be referred to henceforth as level 2 MR personnel. It is the responsibility of the MR medical director not only to identify the necessary training, but also to identify those individuals who qualify as level 2 MR personnel. It is understood that the medical director will have the necessary education and experience in MR safety to qualify as level 2 MR personnel. (See Appendix 1.)
- c. All those not having successfully complied with this MR safety instruction guideline shall be referred to henceforth as non-MR personnel. Specifically, non-MR personnel will be the terminology used to refer to any individual or group who has not within the previous 12 months undergone the designated formal training in MR safety issues defined by the MR safety director of that installation.

3. Patient and non-MR personnel screening

- a. All non-MR personnel wishing to enter Zone III must first pass an MR safety screening process. Only MR personnel are authorized to perform an MR safety screen before permitting non-MR personnel into Zone III.
- b. The screening process and screening forms for patients, non-MR personnel, and MR personnel should be essentially identical. Specifically, one should assume that non-MR personnel, health care practitioners, or MR personnel may enter the bore of the MR imager during the MR imaging process.

Examples of this might include when a pediatric patient cries for his mother, who then leans into the bore, or when the anesthetist leans into the bore to manually ventilate a patient in the event of a problem.

c. Metal detectors

The usage in MR environments of conventional metal detectors which do not differentiate between ferrous and nonferromagnetic materials is not recommended. Reasons for this recommendation against conventional metal detector usage include, among others:

1. They have varied—and variable—sensitivity settings.
2. The skills of the operators can vary.
3. Today's conventional metal detectors cannot detect, for example, a 2 × 3 mm, potentially dangerous ferromagnetic metal fragment in the orbit or near the spinal cord or heart.
4. Today's conventional metal detectors do not differentiate between ferromagnetic and nonferromagnetic metallic objects, implants, or foreign bodies.
5. Metal detectors should not be necessary for the detection of large metallic objects, such as oxygen tanks on the gurney with the patients. These objects are fully expected to be detected—and physically excluded—during the routine patient screening process.

However, ferromagnetic detection systems are currently available that are simple to operate, capable of detecting even very small ferromagnetic objects external to the patient, and now, for the first time, differentiating between ferromagnetic and nonferromagnetic materials. While the use of conventional metal detectors is not recommended, the use of **ferromagnetic detection systems** is recommended as an adjunct to thorough and conscientious screening of persons and devices approaching Zone IV. It should be reiterated that their use is in no way meant to replace a thorough screening practice, which rather should be supplemented by their usage.

- d. Non-MR personnel should be accompanied by, or under the immediate supervision of and in visual or verbal contact with, one specifically identified level 2 MR person for the entirety of their duration within Zone III or Zone IV restricted regions. However, it is acceptable to have non-MR personnel in a changing room or restroom in Zone III without visual contact as long as the personnel and the patient can communicate verbally with each other.

Level 1 MR personnel are permitted unaccompanied access throughout Zones III and IV. Level 1 MR personnel are also explicitly permitted to be responsible for accompanying non-MR personnel into and throughout Zone III, excluding Zone IV. However, level 1 MR personnel are *not* permitted to directly admit, or be designated responsible for, non-MR personnel in Zone IV.

In the event of a shift change, lunch break, etc., no level 2 MR personnel shall relinquish their responsibility to supervise non-MR personnel still within Zone III or Zone IV until such supervision has been formally transferred to another of the site's level 2 MR personnel.

- e. Nonemergent patients should be MR safety-screened on site by a minimum of 2 separate individuals. At least one of these individuals should be level 2 MR personnel. At least one of these 2 screenings should be performed verbally or interactively.

Emergent patients and their accompanying non-MR personnel may be screened only once, providing the screening individual is level 2 MR personnel. There should be no exceptions to this.

- f. Any individual undergoing an MR procedure must remove all readily removable metallic personal belongings and devices on or in them (e.g., watches, jewelry, pagers, cell phones, body piercings [if removable], contraceptive diaphragms, metallic

Safe MR Practices

drug delivery patches [see section I, below], cosmetics containing metallic particles [such as eye make-up], and clothing items that may contain metallic fasteners, hooks, zippers, loose metallic components, or metallic threads). It is therefore advisable to require that the patients or research subjects wear a site-supplied gown with no metal fasteners when feasible.

- g. All patients and non-MR personnel with a history of potential ferromagnetic foreign object penetration must undergo further investigation prior to being permitted entrance to Zone III. Examples of acceptable methods of screening include patient history, plain X-ray films, prior CT or MR studies of the questioned anatomic area, or access to written documentation as to the type of implant or foreign object that might be present. Once positive identification has been made as to the type of implant or foreign object that is within a patient, best-effort assessments should be made to identify the MR compatibility or MR safety of the implant or object. Efforts at identification might include written records of the results of formal testing of the implant prior to implantation (preferred), product labeling regarding the implant or object, and review of peer-reviewed publications regarding MR compatibility and MR safety testing of the make, model, and type of the object. MR safety testing would be of value only if the object or device had not been altered since such testing results had been published.

All patients who have a history of orbit trauma by a potential ferromagnetic foreign body *for which they sought medical attention* are to have their orbits cleared either by plain X-ray or bit films (2 views) [12, 13] or by a radiologist's review and assessment of contiguous cut prior CT or MR images (obtained since the suspected traumatic event), if available.

- h. Conscious, nonemergent patients and research and volunteer subjects are to complete written MR safety screening questionnaires prior to their introduction to Zone III. Family or guardians of nonresponsive patients or of patients who cannot reliably provide their own medical histories are to complete a written MR safety screening questionnaire prior to their introduction to Zone III. These completed questionnaires are then to be reviewed orally with the patient, guardian, or research subject in their entirety prior to permitting the patient or research subject to be cleared into Zone III.

The patient, guardian, or research subject as well as the screening MR staff member must both sign the completed form. This form should then become part of the patient's medical record. No empty responses will be accepted—each question *must* be answered with a “yes” or “no” or specific further information must be provided as requested. A sample pre-MR screening form is provided (see Appendix 3). This is the minimum information to be obtained; more may be added if the site so desires.

- i. Screening of the patient or non-MR personnel with, or suspected of having, an intracranial aneurysm clip should be performed as per the separate MR safe practice guideline addressing this particular topic (see section M, below).
- j. Screening of patients for whom an MR examination is deemed clinically indicated or necessary, but who are unconscious or unresponsive, who cannot provide their own reliable histories regarding prior possible exposures to surgery, trauma, or metallic foreign objects, and for whom such histories cannot be reliably obtained from others:

1. If no reliable patient metal exposure history can be obtained, and if the requested MR examination cannot reasonably wait until a reliable history might be obtained, it is recommended that such patients be physically examined by level 2 MR personnel. All areas of scars or deformities that might be anatomically indicative of an implant, such as on the chest or spine region, and whose origins are unknown and which may have been caused by ferromagnetic foreign bodies, implants, etc., should be subject to plain-film radiography (if recently obtained plain films or CT or MR studies of such areas are not already available). The investigation described above should be made to ensure there are no potentially harmful embedded or implanted metallic foreign objects or devices. All such patients should also undergo plain film imaging of the skull or orbits and chest to exclude metallic foreign objects (if recently obtained plain films or CT or MR studies of such areas are not already available).
2. Monitoring of patients in the MR scanner is sometimes necessary. The potential for thermal injury from excessive RF power deposition exists. Sedated, anesthetized, or unconscious patients may not be able to express symptoms of such injury. This potential for injury is greater on especially higher-field whole-body scanners (e.g., 1 Tesla and above). Distortion of the electrocardiogram within the magnetic field makes interpretation of the ECG complex unreliable, even with filtering used by contemporary monitoring systems. However, routine monitoring of heart rate and rhythm may be accomplished using pulse oximetry, which also eliminates the risks of thermal injury from electrocardiography. Patients who require ECG monitoring and who are unconscious, sedated, or anesthetized should be examined after each imaging sequence, with potential repositioning of the ECG leads and any other electrically conductive material with which the patient is in contact. Alternatively, cold compresses or ice packs could be placed upon all necessary electrically conductive material that touches the patient during scanning.
- k. Final determination of whether or not to scan any given patient with any given implant, foreign body, etc., is to be made by the level 2 MR personnel—designated attending MR radiologist, the MR medical director, or specifically designated level 2 MR personnel following criteria for acceptability predetermined by the medical director.

For implants that are strongly ferromagnetic, an obvious concern is that of magnetic translational and rotational forces upon the implant which might move or dislodge the device from its implanted position. If an implant has demonstrated weak ferromagnetic forces on formal testing, it might be prudent to wait several weeks for fibrous scarring to set in, as this may help anchor the implant in position and help it resist such weakly attractive magnetic forces that might arise in MR environments.

For all implants that have been demonstrated to be nonferrous in nature, however, the risk of implant motion is essentially reduced to those resulting from Lenz's forces alone. These tend to be quite trivial for typical metallic implant sizes of a few centimeters or less. Thus, a waiting period for fibrous scarring to set in is far less important, and the advisability for such a waiting period may well be easily outweighed by the potential clinical benefits of undergoing an MR examination at that time. As always, clinical assessment of the risk–benefit ratio for the par-

ticular clinical situation and patient at hand are paramount for appropriate medical decision making in these scenarios.

It is possible that during the course of an MRI examination an unanticipated ferromagnetic implant or foreign body is discovered within a patient or research subject undergoing the examination. This is typically suspected or detected by means of a sizable field-distorting artifact seen on spin-echo imaging techniques that grows more obvious on longer TE studies and expands markedly on typical moderate or long TE gradient-echo imaging sequences. In such cases, it is imperative that the medical director, safety officer, and/or physician in charge be immediately notified of the suspected findings. This individual should then assess the situation, review the imaging information obtained, and decide what the best course of action might be.

It should be noted that there are numerous potentially acceptable courses that might be recommended which in turn depend upon many factors, including the status of the patient, the location of the suspected ferromagnetic implant/foreign body relative to local anatomic structures, the mass of the implant, etc. Appropriate courses of action might include proceeding with the scan under way, immobilizing the patient and the immediate removal of the patient from the scanner, or other intermediate steps. Regardless of the course of action selected, it is important to note that the forces on the implant will change, and may actually increase, during the attempt to remove the patient from the scanner bore. Further, the greater the rate of motion of the patient/device through the magnetic fields of the scanner bore, the greater the forces acting upon that device will likely be. Thus, it is prudent to ensure that, if at all possible, immobilization of the device during patient extraction from the bore, and the slow, cautious, deliberate rate of extricating the patient from the bore, will likely result in weaker and potentially less harmful forces on the device as it traverses the various static magnetic field gradients associated with the MR imager.

It is also worthy of note that the magnetic fields associated with the MR scanner are distributed throughout space three-dimensionally. Thus, especially for superconducting systems, one should avoid the temptation to have the patient sit up as soon as he or she is physically out of the bore. Doing so may expose the ferrous object to still-significant torque- and translation-related forces despite the patient's being physically outside the scanner bore. It is therefore advisable to continue to extract the patient along a straight line course parallel to the center of the magnet while the patient remains immobilized until they are as far as physically possible from the MR imager itself, before any other patient/object motion vector is attempted or permitted.

1. All non-MR personnel (e.g., patients, volunteers, varied site employees, and professionals) with implanted cardiac pacemakers, autodefibrillators, diaphragmatic pacemakers, or other electro-mechanically activated devices upon which the non-MR personnel is dependent should be precluded from Zone IV and physically restrained from the 5-gauss line unless specifically cleared in writing by a level 2 MR personnel-designated attending radiologist or the medical director of the MR site. In such circumstances, a specific defending risk-benefit rationale should be provided in writing and signed by the authorizing radiologist.

Should it be determined that non-MR personnel wishing to accompany a patient into an MR scan room require their orbits to be

cleared by plain-film radiography, a radiologist must first discuss with the non-MR personnel that plain X-ray films of their orbits are required prior to permitting them access to the MR scan room. Should they still wish to proceed with access to Zone IV or within the 5-gauss line, and should the attending radiologist deem it medically advisable that they do so (e.g., for the care of their child about to undergo an MR study), written informed consent should be provided by these accompanying non-MR personnel prior to their undergoing X-ray examination of their orbits.

- m. MR scanning of patients, prisoners, or parolees with metallic prisoner-restraining devices or RF ID or tracking bracelets could lead to theoretic adverse events, including: (1) ferromagnetic attractive effects and resultant patient injury, (2) possible ferromagnetic attractive effects and potential damage to the device or its battery pack, (3) RF interference with the MRI study and secondary image artifact, (4) RF interference with the functionality of the device, (5) RF power deposition and heating of the bracelet or tagging device or its circuitry and secondary patient injury (if the bracelet were in the anatomic volume of the RF transmitter coil being used for imaging). Therefore, when requested to scan a patient, prisoner, or parolee wearing RF bracelets or metallic handcuffs or ankle cuffs, request that the patient be accompanied by the appropriate authorities who can and will remove the restraining device prior to the MR study and be charged with its replacement following the examination.
- n. Firefighter, police, and security safety considerations: For the safety of firefighters and other emergent services responding to an emergent call at the MR site, it is recommended that all fire alarms, cardiac arrests, or other emergent service response calls originating from or located in the MR site should be forwarded simultaneously to a specifically designated individual from among the site's MR personnel. This individual should, if possible, be on site prior to the arrival of the firefighters or emergent responders to ensure that they do not have free access to Zone III or Zone IV. The site might consider assigning appropriately trained security personnel, who have been trained and designated as MR personnel, to respond to such calls.

In any case, all MR sites should arrange to prospectively educate their local fire marshals, firefighters' associations, and police or security personnel about the potential hazards of responding to emergencies in the MR suite.

It should be stressed that even in the presence of a true fire (or other emergency) in Zone III or Zone IV, the magnetic fields may be present and fully operational. Therefore, free access to Zone III or Zone IV by firefighters or other non-MR personnel with air tanks, axes, crowbars, other firefighting equipment, guns, etc., might prove catastrophic or even lethal to those responding or to others in the vicinity.

As part of the Zone III and Zone IV restrictions, all MR sites must have clearly marked, readily accessible MR-conditional or MR-safe fire extinguishing equipment physically stored in Zone III or Zone IV. All conventional fire extinguishers and other firefighting equipment not tested and verified safe in the MR environment should be restricted from Zone III.

For superconducting magnets, the helium (and the nitrogen as well, in older MR magnets) is not flammable and does not pose a fire hazard directly. However, the liquid oxygen that can result from the supercooled air in the vicinity of the released

Safe MR Practices

gases might well increase the fire hazard in this area. If there are appropriately trained and knowledgeable MR personnel available during an emergency to ensure that emergency response personnel are kept out of the MR scanner or magnet room and away from the 5-gauss line, quenching the magnet during a response to an emergency or fire should not be a requirement.

However, if the fire is in such a location where Zone III or Zone IV needs to be entered for whatever reason by firefighting or emergency response personnel and their firefighting and emergent equipment, such as air tanks, crowbars, axes, and defibrillators, a decision to quench a superconducting magnet should be *very* seriously considered to protect the health and lives of the emergent responding personnel. Should a quench be performed, appropriately designated MR personnel still need to ensure that *all* non-MR personnel (including and especially emergent response personnel) continue to be restricted from Zones III and IV until the designated MR personnel has personally verified that the static field is either no longer detectable or at least sufficiently attenuated as to no longer present a potential hazard to one moving by it with, for example, large ferromagnetic objects such as air tanks or axes.

For resistive systems, the magnetic field of the MR scanner should be shut down as completely as possible and verified as such prior to permitting the emergency response personnel access to Zone IV. For permanent, resistive, or hybrid systems whose magnetic fields cannot be completely shut down, MR personnel should ideally be available to warn the emergency response personnel that a very powerful magnetic field is still operational in the magnet room.

4. MR personnel screening

All MR personnel are to undergo an MR screening process as part of their employment interview process to ensure their safety in the MR environment. For their own protection and for the protection of the non-MR personnel under their supervision, all MR personnel must immediately report to the MR medical director any trauma, procedure, or surgery they experience or undergo in which a ferromagnetic metallic object or device may have become introduced within or on them. This will permit appropriate screening to be performed on the employee to determine the safety of permitting that employee into Zone III.

5. Device and object screening

Ferrous objects, including those brought by patients, visitors, contractors, etc., should be restricted from entering Zone III, whenever practical.

As part of the Zone III site restriction and equipment testing and clearing responsibilities, all sites should have ready access to a strong handheld magnet (≥ 1000 gauss). This will enable the site to test external, and even some superficial internal, devices or implants for the presence of grossly detectable ferromagnetic attractive forces.

- a. All portable metallic or partially metallic devices that are on or external to the patient (e.g., oxygen cylinders) are to be positively identified in writing as ferromagnetic or, alternatively, nonferromagnetic and safe or conditionally safe in the MR environment prior to permitting them into Zone III. For all device or object screening, verification and positive identification should be in writ-

ing. Examples of devices that need to be positively identified include fire extinguishers, oxygen tanks, and aneurysm clips.

- b. External devices or objects demonstrated to be ferromagnetic and MR unsafe or incompatible in the MR environment may still, under specific circumstances, be brought into Zone III if, for example, they are deemed by MR personnel to be necessary and appropriate for patient care. They should only be brought into Zone III if they are under the direct supervision of specifically designated level 1 or level 2 MR personnel who are thoroughly familiar with the device, its function, and the reason supporting its introduction to Zone III. The safe utilization of these devices while they are present in Zone III will be the responsibility of specifically named level 1 or 2 MR personnel. These devices must be appropriately physically secured or restricted at all times during which they are in Zone III to ensure that they do not inadvertently come too close to the MR scanner and accidentally become exposed to static magnetic fields or gradients that might result in their becoming either hazardous projectiles or no longer accurately functional.
- c. Never assume MR compatibility or safety information about the device if it is not clearly documented in writing. All unknown external objects or devices being considered for introduction beyond Zone II should be tested with a strong handheld magnet (≥ 1000 gauss) for ferromagnetic properties before permitting them entry to Zone III. The results of such testing, as well as the date, time, and name of the tester, and methodology used for that particular device, should be documented in writing. If a device has not been tested, or if its MR compatibility or safety status is unknown, it should *not* be permitted unrestricted access to Zone III.
- d. All portable metallic or partially metallic objects that are to be brought into Zone IV must be properly identified and appropriately labeled utilizing the current FDA labeling criteria developed by ASTM (American Society for Testing and Materials) International (<http://www.astm.org>) (see Fig. 2). Those items which are wholly nonmetallic should be identified with a square green "MR safe" label. Items which are clearly ferromagnetic should be identified as "not MR safe" and labeled appropriately with the corresponding round red label with a slash through it. Objects with an "MR conditional" rating should be affixed with a triangular yellow MR conditional label prior to being taken into the scan room/Zone IV.

As noted in the introduction to this section B.5, above, if MR safety data are not prospectively available for a given device, initial testing for the purpose of this labeling is to be accomplished by the site's MR personnel by exposing the metallic object to a handheld magnet (≥ 1000 gauss). If grossly detectable attractive forces are observed between the object being tested or any of its components and the handheld magnet, it is to be labeled with a circular red "not MR safe" label. If no or negligible attractive forces are observed, a triangular yellow "MR conditional" label is to be attached to the object. It is only when the composition of an object and its components are known to be nonmetallic that the green "MR safe" label is to be affixed to a device or object.

Particularly with regard to nonclinical and incidental equipment, current products marketed with ill-defined terminology such as "non-magnetic," or outdated classifications such as "MR-compatible," should not be presumed to conform to a particular current ASTM classification. Similarly, any product marketed as "MR safe" but with metallic construction or components should be treated with suspicion. Objects intended for

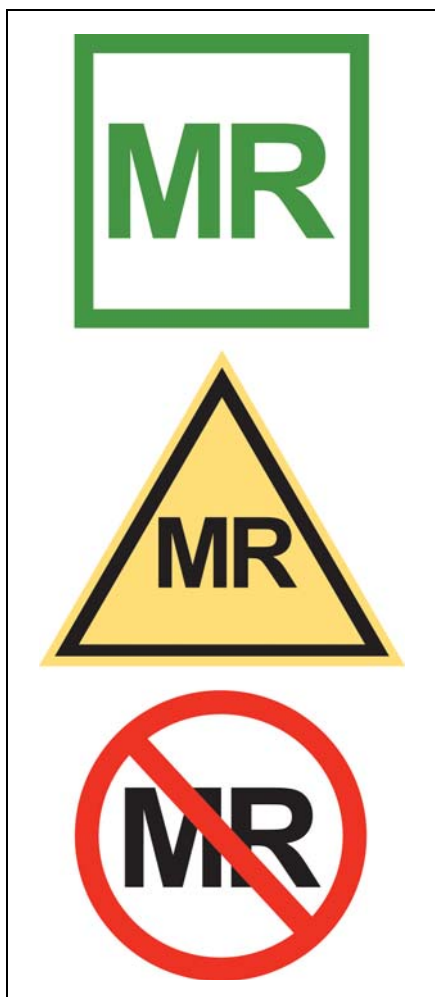


Fig. 2—U.S. Food and Drug Administration labeling criteria (developed by ASTM [American Society for Testing and Materials] International) for portable objects taken into Zone IV. Square green “MR safe” label is for wholly nonmetallic objects, triangular yellow label is for objects with “MR conditional” rating, and round red label is for “not MR safe” objects.

use in Zone IV, including nonclinical incidental products such as stepping stools or ladders, which are not provided with manufacturer or third-party MR safety test results under the new ASTM criteria, should be site tested as described above.

- e. Decisions based on published MR compatibility or safety claims should recognize that all such claims apply to specifically tested static field and static gradient field strengths—for example, “MR conditional, having been tested to be safe up to 3.0 Tesla at gradient strengths of 400 G/cm,” or “MR conditional, having been tested to be safe up to 1.5 Tesla up to maximum static gradient fields experienced in an unshielded 1.5-Tesla [manufacturer’s name] whole-body MR scanner tested 1.5 feet (roughly 45 cm) within the bore.”
- f. It should be noted that alterations performed by the site on MR safe, MR unsafe, and MR conditional equipment or devices may alter the MR safety or compatibility properties of the device. For example, tying a ferromagnetic metallic twisting binder onto a sign labeling the device as MR conditional or MR safe might result in artifact induction—or worse—if introduced into the MR scanner.

Lenz’s Forces:

Faraday’s law states that a moving or changing magnetic field will induce a voltage in a perpendicularly oriented electrical conductor. Lenz’s law builds upon this and states that the induced voltage will itself be such that it will secondarily generate its own magnetic field whose orientation and magnitude will oppose those of the initial time-varying magnetic field that created it in the first place. For example, if an electrical conductor is moved perpendicularly toward the magnetic field, B_0 , of an MR scanner, even if this conductor is not grossly ferromagnetic, the motion itself will result in the generation of voltages in this conductor whose magnitude is directly proportional to the rate of motion as well as the spatial gradient of the magnetic field, B_0 , through which it is being moved. Conducting objects turning in the static field will also experience a torque due to the induced eddy currents. Lenz’s law states that this induced current will in turn create a magnetic field whose orientation will oppose the B_0 magnetic field that created this current.

Thus, moving a large metallic but nonferromagnetic electrical conductor toward the magnet bore will result in the induction of a voltage and associated magnetic field which will orient in such a manner and at such a strength as to oppose the motion of the metallic object into the bore of the MR scanner. If, for example, one tries to move a nonferrous oxygen tank into the bore of an MR scanner, as the scanner bore is approached Lenz’s forces will be sufficiently strong to virtually stop forward progress of the tank. Further, the faster one moves the tank into the bore, the greater the opposing force that is created to stop this motion.

This also has potential consequences for large implanted metallic devices such as certain metallic nonferrous infusion pumps. Although they may not pose a projectile hazard, rapid motion of the patient/implant perpendicular to the magnetic field of the MR imager can be expected to result in forces on the implant that would oppose this motion and may likely be detected by the patient. If the patient were to complain of experiencing forces tugging or pulling on the implant, this might bring the patient or health care personnel to erroneously conclude that there were ferrous components to the device, which might lead to cancellation of the examination. Slowly moving such large metallic devices into and out of the bore is a key factor in decreasing any Lenz’s forces that might be induced and in decreasing the likelihood of a misunderstanding or an unnecessary study cancellation.

C. MR Technologists

1. MR technologists should be ARRT (American Registry of Radiologic Technologists)—registered technologists (RTs). Furthermore, all MR technologists must be trained as level 2 MR personnel during their orientation prior to being permitted free access to Zone III.
2. All MR technologists will maintain current certification in American Heart Association basic life support at the health care provider level.
3. Except for emergent coverage, there will be a minimum of 2 MR technologists or one MR technologist and one other individual with the designation of MR personnel in the immediate Zone II through Zone IV MR environment. For emergent coverage, the MR technologist can scan with no other individuals in their Zone II through Zone IV environment as long as there is in-house, ready emergent coverage by designated department of radiology MR personnel (e.g., radiology house staff or attending radiologist).

D. Pregnancy-Related Issues

1. Health care practitioner pregnancies

Pregnant health care practitioners are permitted to work in and around the MR environment throughout all stages of their pregnancy [14]. Acceptable activities include, but are not limited to, positioning patients, scanning, archiving, injecting contrast material, and entering the MR scan room in response to an emergency. Although permitted to work in and around the MR environment, pregnant health care practitioners are requested not to remain within the MR scanner bore or Zone IV during actual data acquisition or scanning.

2. Patient pregnancies

Present data have not conclusively documented any deleterious effects of MR imaging exposure on the developing fetus. Therefore, no special consideration is recommended for the first, versus any other, trimester in pregnancy. Nevertheless, as with all interventions during pregnancy, it is prudent to screen women of reproductive age for pregnancy prior to permitting them access to MR imaging environments. If pregnancy is established, consideration should be given to reassessing the potential risks versus benefits of the pending study in determining whether performance of the requested MR examination could safely wait until the end of the pregnancy.

- a. Pregnant patients can be accepted to undergo MR scans at any stage of pregnancy if, in the determination of a level 2 MR personnel-designated attending radiologist, the risk-benefit ratio to the patient warrants that the study be performed. The radiologist should confer with the referring physician and document the following in the radiology report or the patient's medical record:
 1. The information requested from the MR study cannot be acquired via nonionizing means (e.g., ultrasonography).
 2. The data are needed to potentially affect the care of the patient or fetus *during* the pregnancy.
 3. The referring physician does not feel it is prudent to wait until the patient is no longer pregnant to obtain these data.
- b. MR contrast agents should *not* be routinely provided to pregnant patients. This decision, too, is one that must be made on a case-by-case basis by the covering level 2 MR personnel-designated attending radiologist who will assess the risk-benefit ratio for that particular patient.

The decision to administer a gadolinium-based MR contrast agent to pregnant patients should be accompanied by a well-documented and thoughtful risk-benefit analysis. This analysis should be able to defend a decision to administer the contrast agent based on overwhelming potential benefit to the patient or fetus outweighing the theoretic but potentially real risks of long-term exposure of the developing fetus to free gadolinium ions.

Studies have demonstrated that gadolinium-based MR contrast agents pass through the placental barrier and enter the fetal circulation. From there, they are filtered in the fetal kidneys and then excreted into the amniotic fluid. In this location the gadolinium-chelate molecules are in a relatively protected space and may remain in this amniotic fluid for an indeterminate amount of time before finally being reabsorbed and eliminated. As with any equilibrium situation involving any dissociation constant, the longer the chelate molecule remains in this space, the greater the potential for disso-

ciation of the potentially toxic gadolinium ion from its chelate molecule. It is unclear what impact such free gadolinium ions might have if they were to be released in any quantity in the amniotic fluid. Certainly, deposition into the developing fetus would raise concerns of possible secondary adverse effects.

The risk to the fetus with administration of gadolinium-based MR contrast agents remains unknown and may be harmful.

- c. It is recommended that pregnant patients undergoing an MR examination provide written informed consent documenting that they understand the potential risks and benefits of the MR procedure to be performed, are aware of the alternative diagnostic options available to them (if any), and wish to proceed.

E. Pediatric MR Safety Concerns

1. Sedation and monitoring issues

Children form the largest group requiring sedation for MRI, largely because of their inability to remain motionless during scans. Sedation protocols may vary from institution to institution according to the procedures performed (diagnostic vs interventional), the complexity of the patient population (healthy preschoolers vs premature infants), the method of sedation (mild sedation vs general anesthesia), and the qualifications of the sedation provider.

Adherence to standards of care mandates following the sedation guidelines developed by the American Academy of Pediatrics [15, 16], the American Society of Anesthesiologists [17], and the Joint Commission on Accreditation of Healthcare Organizations (JCAHO) [18]. In addition, sedation providers must comply with protocols established by the individual state and the practicing institution. These guidelines require the following provisions:

- a. Preprocedural medical history and examination for each patient
- b. Fasting guidelines appropriate for age
- c. Uniform training and credentialing for sedation providers
- d. Intraprocedural and postprocedural monitors with adaptors appropriately sized for children (compatible with the magnetic field)
- e. Method of patient observation (window, camera)
- f. Resuscitation equipment, including oxygen delivery and suction
- g. Uniform system of record keeping and charting (with continuous assessment and recording of vital signs)
- h. Location and protocol for recovery and discharge
- i. Quality assurance program that tracks complications and morbidity.

For the neonatal and the young pediatric population, special attention is needed in monitoring body temperature for both hypo- and hyperthermia in addition to other vital signs [19]. Temperature-monitoring equipment that is approved for use in the MR suite is becoming more readily available. Commercially available, MR-approved neonatal isolation transport units and other warming devices are also available for use during MR scans.

2. Pediatric screening issues

Children may not be reliable historians and, especially in cases of older children and teenagers, should be questioned both in the presence of parents or guardians and separately to maximize the possibility that all potential dangers are disclosed. Therefore, it is recommended that children be gownned before entering Zone IV to help ensure that no metallic objects, toys, etc. inadvertently find

their way into Zone IV. Pillows, stuffed animals, and other comfort items brought from home represent real risks and should be discouraged from entering Zone IV. If unavoidable, each such item should be carefully checked with the powerful handheld magnet and perhaps again in the MR scanner prior to permitting the patient to enter Zone IV with the object in order to ensure that it does not contain any objectionable metallic components.

3. MR safety of accompanying family or personnel

Although any age patient might request that others accompany them for their MR examination, this is far more common in the pediatric population. Those accompanying or remaining with the patient should be screened using the same criteria as anyone else entering Zone IV.

In general, it would be prudent to limit accompanying adults to a single individual. Only a qualified, responsible MR physician should make screening criteria exceptions.

Hearing protection and MR safe/MR conditional seating are recommended for accompanying family members within the MR scan room.

F. Time-Varying Gradient Magnetic Field–Related Issues: Induced Voltages

Types of patients needing extra caution:

Patients with implanted or retained wires in anatomically or functionally sensitive areas (e.g., myocardium or epicardium, implanted electrodes in the brain) should be considered to be at higher risk, especially from faster MRI sequences, such as echo-planar imaging (which may be used in such sequences as diffusion-weighted imaging, functional imaging, perfusion-weighted imaging, MR angiographic imaging, etc.). The decision to limit the dB/dt (rate of magnetic field change) and maximum strength of the magnetic field of the gradient subsystems during imaging of such patients should be reviewed by the level 2 MR personnel—designated attending radiologist supervising the case or patient.

G. Time-Varying Gradient Magnetic Field–Related Issues: Auditory Considerations

1. All patients and volunteers should be offered and encouraged to use hearing protection prior to undergoing any imaging in the MR scanners.
2. All patients or volunteers in whom research sequences are to be performed (i.e., MR scan sequences that have not yet been approved by the Food and Drug Administration) are to have hearing protective devices *in place* prior to the initiation of any MR sequences. Without hearing protection in place, MRI sequences that are not FDA-approved should not be performed on patients or volunteers.

H. Time-Varying Radiofrequency Magnetic Field–Related Issues: Thermal

1. All unnecessary or unused electrically conductive materials should be removed from the MR system before the onset of imaging. It is not sufficient to merely to “unplug” or disconnect unused, unnecessary electrically conductive material and leave it within the MR scanner with the patient during imaging. All electrical connections, such as on surface coil leads or monitoring devices, must be visu-

ally checked by the scanning MR technologist prior to each use to ensure the integrity of the thermal and electrical insulation.

2. Electrical voltages and currents can be induced in electrically conductive materials that are within the bore of the MR imager during the MR imaging process. This might result in the heating of this material by resistive losses. This heat might be of a caliber sufficient to cause injury to human tissue. Among the variables that determine the amount of induced voltage or current is the consideration that the larger the diameter of the conductive loops, the greater the potential induced voltages or currents, and thus the greater the potential for resultant thermal injury to adjacent or contiguous patient tissue.

Therefore, when electrically conductive material (wires, leads, implants, etc.) are required to remain within the bore of the MR scanner with the patient during imaging, care should be taken to ensure that no large-caliber electrically conducting loops (including patient tissue; see section H.5, below) are formed within the MR scanner during imaging. Furthermore, it is possible, with the appropriate configuration, lead length, static magnetic field strength, and other settings, to introduce resonant circuitry between the transmitted RF power and the lead. This could result in very rapid and clinically significant lead heating, especially at the lead tips, in a matter of seconds to a magnitude sufficient to result in tissue thermal injury or burns. This can also theoretically occur with implanted leads or wires, even when they are not connected to any other device at either end. For illustration, the FDA has noted several reports of serious injury, including coma and permanent neurologic impairment, in patients with implanted neurologic stimulators who underwent MR imaging examinations. The injuries in these instances resulted from heating of the electrode tips [20, 21].

Further, it is entirely possible for a lead or wire to demonstrate no significant heating while undergoing MR imaging examinations at 1.5 Tesla, yet demonstrate clinically significant and potentially harmful degrees of heating within seconds at, for example, 3 Tesla. It has also been demonstrated that leads may show no significant heating at 3 Tesla yet may rapidly heat to hazardous levels when undergoing MR imaging at, for example, 1.5 Tesla (personal observation, MR safety testing, E. Kanal, MD, University of Pittsburgh Medical Center MR Research Center, 8/28/05). Thus, at no time should a label of “MR conditionally safe for thermal issues at [a given field strength]” be applied to any field strength, higher or lower, other than the specific one at which safety was demonstrated.

Thus, exposure of electrically conductive leads or wires to the RF transmitted power during MR scanning should only be performed with caution and with appropriate steps taken to ensure significant lead or tissue heating does not result (see section H.9, below).

3. When electrically conductive materials are required to be within the bore of the MR scanner with the patient during imaging, care should be taken to place thermal insulation (including air, pads, etc.) between the patient and the electrically conductive material, while simultaneously attempting (as much as feasible) to keep the electrical conductor from directly contacting the patient during imaging. It is also appropriate to try to position the leads or wires as far as possible from the inner walls of the MR scanner if the body coil is being used for RF transmission. When it is necessary that electrically conductive leads directly contact the patient during imaging, consideration should be given to prophylactic application of cold compresses or ice packs to such areas.
4. Depending on specific magnet designs, care may be needed to ensure that the patient’s tissue(s) do not directly come into contact with the inner bore of the MR imager during the MRI process. This

is especially important for several higher-field MR scanners. The manufacturers of these devices provide pads and other such insulating devices for this purpose, and manufacturer's guidelines should be strictly adhered to for these units.

5. It is important to ensure the patient's tissues do not form large conductive loops. Therefore, care should be taken to ensure that the patient's arms or legs are not positioned in such a way as to form a large-caliber loop within the bore of the MR imager during the imaging process. For this reason, it is preferable that patients be instructed not to cross their arms or legs in the MR scanner. We are also aware of unpublished reports of thermal injuries that seem to have been associated with skin folds, such as in the region of the inner thighs. While the cause of this is not yet fully understood, it might be prudent to consider ensuring that skin folds and other such examples of tissue-to-tissue contact are minimized or eliminated in the region undergoing radiofrequency energy irradiation.
6. Skin staples and superficial metallic sutures: Patients requested to undergo MR studies in whom there are skin staples or superficial metallic sutures (SMS) may be permitted to undergo the MR examination if the skin staples or SMS are not ferromagnetic and are not in the anatomic volume of RF power deposition for the study to be performed. If the nonferromagnetic skin staples or SMS are within the volume to be RF-irradiated for the requested MR study, several precautions are recommended.
 - a. Warn the patient and make sure that they are especially aware of the possibility that they may experience warmth or even burning along the skin staple or SMS distribution. The patient should be instructed to report immediately if they experience warmth or burning sensations during the study (and not, for example, wait until the "end of the knocking noise").
 - b. It is recommended that a cold compress or ice pack be placed along the skin staples or SMS if this can be safely clinically accomplished during the MRI examination. This will help to serve as a heat sink for any focal power deposition that may occur, thus decreasing the likelihood of a clinically significant thermal injury or burn to adjacent tissue.
7. For patients with extensive or dark tattoos, including tattooed eyeliner, in order to decrease the potential for RF heating of the tattooed tissue, it is recommended that cold compresses or ice packs be placed on the tattooed areas and kept in place throughout the MRI process if these tattoos are in the volume in which the body coil is being used for RF transmission. This approach is especially appropriate if fast spin-echo (or other high RF duty cycle) MRI sequences are anticipated in the study. If another coil is being used for RF transmission, a decision must be made if high RF transmitted power is to be anticipated by the study protocol design. If so, then the above precautions should be followed. Additionally, patients with tattoos that had been placed within 48 hours prior to the pending MR examination should be advised of the potential for smearing or smudging of the edges of the freshly placed tattoo.
8. In the unconscious or unresponsive patient, all attached leads that will be in or partly in the volume undergoing RF irradiation should be covered with a cold compress or ice pack at the lead attachment site for the duration of the MR study.
9. As noted above, it has been demonstrated that resonant circuitry can be established during MRI between the RF energies being transmitted and specific lengths of long electrically conductive wires or leads, which can thus act as efficient antennae. This can result in heating of the tips of these wires or leads to temperatures in excess of 90°C in a

few seconds. Therefore, patients in whom there are long electrically conductive leads, such as Swan-Ganz thermodilation cardiac output-capable catheters or Foley catheters with electrically conductive leads, should be considered at risk for MR studies if the body coil is to be used for RF transmission over the region of the electrically conductive lead. This is especially true for higher-field systems and for imaging protocols utilizing fast spin-echo or other high RF duty cycle MRI sequences. Each such patient should be reviewed and cleared by an attending level 2 radiologist and a risk-benefit ratio assessment performed prior to permitting them access to the MR scanner.

10. The potential to establish substantial heating is itself dependent on multiple factors, including, among others, the static magnetic field strength of the MR scanner (as this determines the transmitted radiofrequencies [RF] at which the device operates) and the length, orientation, and inductance of the electrical conductor in the RF-irradiated volume being studied. *Virtually any lead lengths can produce substantial heating.* Innumerable factors can affect the potential for tissue heating for any given lead. It is therefore critical to recognize that of all electrically conductive implants, it is specifically wires, or leads, that pose the greatest potential hazard for establishing substantial power deposition/heating considerations.

Another important consideration is that as a direct result of the above, it has already been demonstrated in vitro that heating of certain implants or wires may be clinically insignificant at, for example, 1.5 Tesla but quite significant at 3.0 Tesla. However, it has also been demonstrated that specific implants might show no significant thermal issues or heating at 3.0 Tesla, *but may heat to clinically significant or very significant levels in seconds at, for example, 1.5 Tesla.* Thus, it is important to follow established product MR safety guidelines carefully and precisely, applying them to, and only to, the static magnetic field strengths at which they had been tested. MR scanning at either stronger and/or weaker magnetic field strengths than those tested may result in significant heating where none had been observed at the tested field strength(s).

I. Drug Delivery Patches and Pads

Some drug delivery patches contain metallic foil. Scanning the region of the metallic foil may result in thermal injury [22]. Since removal or repositioning can result in altering of patient dose, consultation with the patient's prescribing physician would be indicated in assessing how to best manage the patient. If the metallic foil of the patch delivery system is positioned on the patient so that it is in the volume of excitation of the transmitting RF coil, the case should be specifically reviewed with the radiologist or physician covering the patient. Alternative options may include placing an ice pack directly on the patch. This solution may still substantially alter the rate of delivery or absorption of the medication to the patient (and be less comfortable to the patient, as well). This ramification should therefore not be treated lightly, and a decision to proceed in this manner should be made by a knowledgeable radiologist attending the patient and with the concurrence of the referring physician as well.

If the patch is removed, a specific staff member should be given responsibility for ensuring that it is replaced or repositioned.

J. Cryogen-Related Issues

1. For superconducting systems, in the event of a system quench, it is imperative that all personnel and patients be evacuated from the MR

scan room as quickly as safely feasible and that the site access be immediately restricted to all individuals until the arrival of MR equipment service personnel. This is especially so if cryogenic gases are observed to have vented partially or completely into the scan room, as evidenced in part by the sudden appearance of white “clouds” or “fog” around or above the MR scanner. As noted in section B.3.n above, it is especially important to ensure that all police and fire response personnel are restricted from entering the MR scan room with their equipment (axes, air tanks, guns, etc.) until it can be confirmed that the magnetic field has been successfully dissipated, because there may still be a considerable static magnetic field present despite a quench or partial quench of the magnet [23].

2. It should be pointed out that room oxygen monitoring was discussed by the MR Blue Ribbon Panel and rejected at this time because the present oxygen monitoring technology was considered by industry experts not to be sufficiently reliable to allow continued operation during situations of power outages, etc.

K. Claustrophobia, Anxiety, Sedation, Analgesia, and Anesthesia

Adult and pediatric patient anxiolysis, sedation, analgesia, and anesthesia for any reason should follow established ACR [24, 25], American Society of Anesthesiologists (ASA) [26–29], and JCAHO standards [29].

L. Contrast Agent Safety

I. Contrast agent administration issues

No patient is to be administered prescription MR contrast agents without orders from a duly licensed physician. Intravenous injection-qualified MR technologists may start and attend to peripheral IV access lines if they have undergone the requisite site-specified training in peripheral IV access and have demonstrated and documented appropriate proficiency in this area. IV-qualified MR technologists may administer FDA-approved gadolinium-based MR contrast agents via peripheral IV routes as a bolus or as a slow or continuous injection as directed by the orders of a duly licensed site physician.

Administration of these agents is to be performed according to the ACR policy. The ACR approves of the injection of contrast material and diagnostic levels of radiopharmaceuticals by certified and/or licensed radiologic technologists and radiologic nurses under the direction of a radiologist or his or her physician designee who is personally and immediately available, if the practice is in compliance with institutional and state regulations. There must also be prior written approval by the medical director of the radiology department or service of such individuals. Such approval process must follow established policies and procedures, and the radiologic technologists and nurses who have been so approved must maintain documentation of continuing medical education related to materials injected and to the procedures being performed [30].

2. Prior contrast agent reaction issues

- a. According to the *ACR Manual on Contrast Media* [31], adverse events after intravenous injection of gadolinium seem to be more common in patients who had previous reactions to an MR contrast

agent. In one study, 16 (21%) of 75 patients who had previous adverse reactions to MR contrast agents reacted to subsequent injections of gadolinium [31]. Patients with asthma also seem to be more likely to have an adverse reaction to the administration of a gadolinium-based MR contrast agent. Patients with allergies also seemed to be at increased risk (~2.0–3.7 times, compared with patients without allergies). Patients who have had adverse reactions to iodinated contrast media are more than twice as likely to have an adverse reaction to gadolinium (6.3% of 857 patients) [31].

- b. At present, there are no well-defined policies for patients who are considered to be at increased risk for having an adverse reaction to MR contrast agents. However, the following recommendations are suggested: Patients who have previously reacted to one MR contrast agent can be injected with another agent if they are restudied, and at-risk patients can be premedicated with corticosteroids and, occasionally, antihistamines [31].
- c. All patients with asthma, a history of allergic respiratory disorders, prior iodinated or gadolinium-based contrast reactions, etc., should be followed more closely as they are at a demonstrably higher risk of adverse reaction.

3. Renal disease, gadolinium-based MR contrast agents, and nephrogenic systemic fibrosis (NSF)

a. Overview:

It has been recently noted that over a 4-year period, 20 patients in Denmark and five in Austria developed a very rare disease that is seen only in patients with severely impaired renal function [32, 33]. Each of these patients had been administered Omniscan (gadodiamide, GE Healthcare), a gadolinium-based MR contrast agent (GBMCA), for an MR imaging or angiographic examination within a few weeks or months prior to the onset of the disease. Roughly 17,500 patients are examined using Omniscan in Denmark each year. Since January 2002, about 400 patients with severely impaired renal function had been examined, of which 20, or 5%, to whom Omniscan had been administered, eventually were diagnosed with this disease in that country.

The disease in question, originally known as nephrogenic fibrosing dermopathy (NFD) and now more widely recognized as nephrogenic systemic fibrosis (NSF), was only first observed in 1997 and formally described in 2000 [34]. It is associated with increased tissue deposition of collagen, often resulting in thickening and tightening of the skin (usually involving predominantly the distal extremities but occasionally also the trunk) and fibrosis that may involve other parts of the body, including the diaphragm, heart, lung, pulmonary vasculature, and skeletal muscles. There is no definitive cure, although some anecdotal reports exist of at least partial responses to various therapies such as plasmapheresis, extracorporeal photopheresis, and thalidomide. There are some data that suggest slowing or even reversal of the disease symptoms may accompany improvements in renal function (especially transplantation). The disease is progressive and can be fulminant in approximately 5% of cases and can even be associated with a fatal outcome. It is generally seen in middle-aged patients but has also been seen in the elderly as well as the pediatric population [35, 36]. There may be a special predilection for patients with concurrent hepatic disease, but this is not yet clearly established [37, 38].

A central registry for NSF patients is maintained at Yale University by Dr. Shawn Cowper, one of the physicians who originally described this disease [39]. At the time of this writing (1/25/07), virtually all reg-

istry cases in which records can be located have at least one known exposure to gadolinium within a few days to months prior to the development of clinical symptoms [37, verbal communication with Dr. Cowper, December 2006].

b. The association with gadolinium-based MR contrast agents (GBMCAs):

Besides the initial reports noted above, in August 2, 2006, researchers from the Copenhagen University Hospital in Denmark published in the *Journal of the American Society of Nephrology* [40] the results of their review of all 13 confirmed cases of NSF, in which they found that all 13 had received Omniscan 2–75 days (median, 25 days) prior to the development of NSF. To quote from their manuscript, “No other exposure/event than gadodiamide that was common to more than a minority of the patients could be identified. These findings indicate that gadodiamide plays a causative role in nephrogenic systemic fibrosis.”

In that article, they also reported that these 13 patients with confirmed NSF were among roughly 370 severe renal disease patients whom they had tracked who had undergone gadodiamide exposure/administration for an MR examination, whereas none of > 430 patients with severe renal disease who had not received a GBMCA developed NSF.

Although this association was initially reported between Omniscan and NSF, there are now multiple submitted MedWatch cases [11] that allege that diagnoses of NSF followed intravenous administration of Magnevist (gadopentetate dimeglumine, Schering) as well as intravenous administration of OptiMARK (gadoversetamide, Mallinckrodt), which are other chelates of GBMCAs. It is clear that the vast majority ($\approx 90\%$) of known cases at this time seem to be clearly associated with Omniscan to a degree that is out of proportion to the relative market shares for these agents [41, 42]. As of January 17, 2007, of the > 100 cases of NSF reported to the FDA MedWatch reporting system, 85 are Omniscan-associated, 21 are Magnevist-associated, six are OptiMARK-associated, none are associated with ProHance (gadoteridol, Bracco Diagnostics), and one is associated with MultiHance (gadobenate dimeglumine, Bracco Diagnostics) (although this same patient also received Omniscan 5 days after their MultiHance MR examination, and subsequently developed NSF) (personal communication, Dr. Melanie Blank, FDA, January 18, 2007). It is also important to recognize the substantial lack of scientific process inherent in the MedWatch reporting system, whose self-reported data can be used at best as general-trend-indicating and typically not for more specific analyses. Nevertheless, the data support the FDA’s concern that this association may exist for the administration of other, or perhaps any of the other, FDA-approved GBMCAs and the subsequent development of NSF. Although there is evidence associating the development of NSF in renal failure patients with only some, but not all, of the FDA-approved GBMCAs to date, prudence dictates that at this time similar concerns be applied to all GBMCAs in this regard until more definitive information is forthcoming on this issue.

c. Causation?

There is a conjecture that suggests that if a causative relationship exists, it may be secondary to accumulation of the gadolinium chelate or free gadolinium in the dependent subcutaneous tissues of the lower and upper extremities (where the disease seems to most often initially manifest itself). Further, if there is free gadolinium

released in any quantity, studies have suggested that it may accumulate in and bind to bone [43]. Very recent initial reports have apparently demonstrated the presence of gadolinium in the biopsies of tissues of NSF patients [44, 45]. In one control individual without NSF, no gadolinium was found using the same electron dispersion spectroscopy technique.

It should also be added that the very visualization of gadolinium in the scanning electron micrographs (SEM) noted in these two recent publications [44, 45] itself is strong evidence that dissociation of the gadolinium from its chelate has occurred. This can be related to the observation that in the process of preparing the tissue for SEM, water-soluble forms of gadolinium would have likely been removed from the specimen, leaving only the insoluble forms to precipitate out (verbal communication, Michael Tweedle, Bracco Diagnostics, January 2007, and Hanns-Joachim Weinmann, Bayer Schering Pharma, January 19, 2007). These precipitates are likely to be largely gadolinium phosphates (verbal communication, Hanns-Joachim Weinmann, January 19, 2007), but this is neither definite nor universally established.

Additionally, it has been noted by several investigators that the development of NSF seemed to most commonly (although not exclusively) follow high-dose administration of gadolinium-based MR contrast agents. This dose–response observation also suggests a possible etiologic role of these agents in the development of NSF in these patients [37].

Although a definitive causal relationship between GBMCA administration to severe renal disease patients and the development of NSF has not been absolutely established, it certainly does appear that gadolinium administration is quite likely a necessary factor in the development of NSF at this time. If a causative role is postulated or even demonstrated, it is unclear whether the causative agent is released free gadolinium, prolonged exposure to abnormally high doses of the gadolinium-plus-chelate molecule, the chelate itself, or some combination of the above with other factors that might be relatively unique to the biochemical milieu of the patient with severe renal failure. This is supported in part by the observation that in several of the publications, including the initial report from the Danish Medicines Agency [33, 37], the incidence of developing NSF in patients with severe or end-stage renal disease after being administered Omniscan appears to be roughly only 3–5%.

There are early data that suggest that elevated levels of phosphate, iron, zinc, or copper [46] or the presence of Fosrenol (lanthanum carbonate, Shire) might serve as efficient competitors for the “attention” of the chelate molecule, so to speak, and increase the concentration of free gadolinium (Gd^{3+}) in the patient, which might therefore increase the potential of the patient to develop NSF. A history of multiple prior GBMCA administrations also seems to be associated with an increased incidence of subsequent development of NSF.

d. Gadolinium toxicity:

Free gadolinium ion exists most commonly in a 3^+ charged form that inhibits those chemical processes that depend upon the influx of calcium (2^+) ions, such as cardiac and skeletal muscle, neurologic discharge, normal coagulation pathways, some enzymatic reactions, etc.

e. FDA guidance:

On December 22, 2006, the FDA issued an update [47] to their earlier (June 9, 2006) public health advisory (PHA) [48]. This new

version has significantly changed from the prior one in several areas. One of these modifications includes the fact that the new version now includes wording that recommends caution in administering GBMCAs to patients with moderate to end-stage renal disease as well as consideration of providing hemodialysis treatment immediately after administration of this agent for patients in this category of renal disease who receive these agents. (The prior advisory recommended caution, especially in patients with end-stage renal disease, with glomerular filtration rates of $< 15 \text{ mL/min/1.73 m}^2$ [48].) Quoting from this more recent advisory [47]:

When a patient with moderate to end-stage kidney disease needs an imaging study, select imaging methods other than MRI or MRA with a gadolinium-based contrast agent for the study whenever possible. If these patients must receive a gadolinium-based contrast agent, prompt dialysis following the MRI or MRA should be considered.

Average excretory rates of gadolinium were 78%, 96%, and 99% in the first to third hemodialysis sessions, respectively, in end-stage renal disease patients who received Magnevist [49]. One study has found that the mean half-life of gadodiamide is 1.3 hours in healthy volunteers, 34.3 hours in patients with a glomerular filtration rate (GFR) range of 2–10 mL/min/1.73 m², 2.6 hours in hemodialysis patients, and 52.7 hours in peritoneal dialysis patients [50]. It is also known that different hemodialysis membranes have been demonstrated to vary in their effectiveness at clearing the administered GBMCA [51].

It should be pointed out that virtually all present data seem to indicate that the vast majority of NSF patients to date had either severe or end-stage renal disease at the time of diagnosis or administration of the GBMCA, with most already being on dialysis. The official National Kidney Foundation staging system classifies patients with glomerular filtration rates between 30 and 59 mL/min/1.73 m² as having stage 3 or moderate chronic kidney disease (CKD), between 15 and 29 mL/min/1.73 m² as stage 4 or severe CKD, and those with GFR $< 15 \text{ mL/min/1.73 m}^2$ or on dialysis as having stage 5 or end-stage CKD. More than one of four adults over age 70 has a GFR of $< 60 \text{ mL/min/1.73 m}^2$, and roughly 7.7 million Americans have a GFR between 30 and 60 mL/min/1.73 m² [52]. Based on NHANES III 1988–1994 (the Third National Health and Nutrition Examination Survey of the CDC) [53], the prevalence of a GFR $< 60 \text{ mL/min/1.73 m}^2$ in US adults ≥ 20 years of age was 8.0%, or more than one of every 13 adults. By age 70, the normal mean value is approximately 70 mL/min/1.73 m². For adults age 70 and older, the prevalence of GFR $< 60 \text{ mL/min/1.73 m}^2$ is roughly 26%, or more than one in four. Finally, the normal GFR for neonates ≤ 8 weeks of age is $< 65 \text{ mL/min/1.73 m}^2$ [54]. Therefore, an advisory statement worded in this manner might result in exposing patients to the potentially greater risks of hemodialysis or in withholding contrast enhancement for their studies. Since the elderly population are among the greatest users of MRI today, this advisory is especially of concern.

f. Other guidance resources:

An overview of this disease, as well as our recommendations for guidelines regarding NSF, renal disease patients, and gadolinium-based MR contrast agent administration, was accepted for publication in *Radiology* and is available for download from *Radiology's* online site [55].

The European Agency for the Evaluation of Medicinal Products (EMA) has recently issued a recommendation [56] to consider the administration of Omniscan (and OptiMARK, although the latter is not licensed in Europe) as contraindicated in patients with severe renal disease (GFR $< 30 \text{ mL/min/1.73 m}^2$) or those who have had or will be undergoing a liver transplant. They also warn that for children up to 1 year of age, because their kidneys are immature, one should be most cautious about administering Omniscan (or OptiMARK). For the other non-Omniscan gadolinium-based MR contrast agents (GBMCAs), they advise simply that there is a possibility of NSF resulting with some GBMCAs in patients with severe renal disease. The European Pharmacovigilance Working Party (PhVWP) and the United Kingdom Commission on Human Medicines (CHM) do not recommend dialysis after administration of GBMCAs, even for hemodialysis patients [56].

As noted above, the FDA continues to recommend considering immediate hemodialysis for any patient with moderate, severe, or end-stage renal disease receiving any GBMCA [47].

g. Recommendations:

At this stage, the following guidelines are recommended when considering administering a GBMCA to patients with renal failure/disease:

The development of NSF in patients with renal disease has followed the administration of some, but not all, of the FDA-approved GBMCAs. To date, the development of NSF has been associated with the isolated prior administration of—especially, and clearly predominantly—Omniscan (at rates that exceed those associated with simple market share), but also Magnevist and OptiMARK. Nevertheless, it is thought to be appropriate to assume for now that a potential association might exist for all five FDA-approved gadolinium-based MR contrast agents until there are more definitive data to suspect otherwise.

At this time, no special treatment or handling is recommended for kidney disease patients with stage 1 or 2 chronic kidney disease (defined as presence of kidney damage with GFR $> 90 \text{ mL/min/1.73 m}^2$ or GFR between 60 and 89 mL/min/1.73 m², respectively). The only exception to this is that patients with any level of renal disease should not receive Omniscan for their contrast-enhanced MR examinations. This is an opinion shared by others [57] and seems prudent for all renal disease patients.

Prospectively checking patient renal function, serum creatinine level, or glomerular filtration rate prior to accepting a patient for an MR imaging or angiographic examination is specifically not required. Among the reasons for this is that roughly 90% of NSF patients seem to already be on dialysis and the majority of the remainder seem to be stage 5 or stage 4. Add to this the fact that one could avoid administering any of the agents with which NSF has been most strongly associated, and the fact that even in patients with severe or end-stage renal disease the incidence of developing NSF seems to be around 3–5%. Therefore, specific prospective hematologic screening is not felt to be warranted. Instead, it is recommended that all requests for MR be pre-screened, with an additional question inquiring about the presence of a history of “kidney disease or dialysis.” If the disease is present but quite mild (stages 1 or 2), modification of how the study should be performed (relative to a patient with no renal disease) does not appear to be indicated except for the avoidance of Omniscan. Conversely, if the disease is present and severe or end-stage in nature, the patient will often be aware of this level of kidney disease and will likely be under

physician care for this condition. *The American Journal of Kidney Diseases* states [54]: "In general, patients with GFR <30 mL/min/1.73 m² should be referred to a nephrologist." Thus, selecting patients with a GFR threshold of roughly 30 mL/min/1.73 m² or already on dialysis (i.e., stages 4 and/or 5) as the level for which special consideration (including possibly hemodialysis) should be given, might represent a medically reasonable approach to, and compromise on, this issue. For patients with stage 3 CKD, the potential risks associated with withholding an MR imaging or angiographic examination could outweigh the potential risk of developing NSF, given the very few number of patients with putative GFR < 60 mL/min/1.73 m² who have been reported to have developed NSF. Further data are clearly needed to clarify the potential risk for stage 3 CKD patients given the few cases reported and the large number of patients with stage 3 CKD and who are predominantly older than age 70 who would be affected.

For all patients with stage 3, 4, or 5 kidney disease or those with acute kidney injury (AKI), it is recommended that one consider refraining from administering any GBMCAs unless a risk-benefit assessment for that particular patient indicates that the benefit of doing so clearly outweighs the potential risk(s). Similar reasoning applies equally to patients with protected regions which the gadolinium chelate might enter but from which it might not be readily cleared. An example of such a space is the amniotic fluid, in which these contrast agents can accumulate shortly after intravenous administration (personal observation and verbal communication, Emanuel Kanal, 1988).

When risk-benefit assessments warrant administration of a GBMCA to patients with stages 3–5 renal disease (moderate to end-stage) or AKI, consideration should be given to administering the lowest dose that would provide the diagnostic benefit being sought, with a half-dose, if clinically acceptable, being considered the default standard dose for such patients. The study should be monitored during its execution and prior to contrast administration to ensure that the administration of the GBMCA is still deemed necessary and indicated at that time. Postponing the examination in patients with AKI until renal function has recovered should also be considered if clinically feasible.

Standard medical practice dictates that for all patients who receive a contrast agent, the type, dose, and route of administration are to be documented in a physician order and in the report. However, patients with moderate to end-stage (stages 3–5) renal disease who are to undergo contrast-enhanced MR imaging examinations of any kind must have a written order to this effect for this agent from the radiologist approving the examination. This order must arise explicitly from the radiologist and NOT from either a referring physician or an MR imaging protocol standing order. The name of the patient, the name and specific brand of GBMCA, dose, route, and rate of administration should all be explicitly specified on the order, along with the date and signature of the requesting radiologist.

Prospective documentation of a risk-benefit assessment for each such patient is considered advisable. It is recommended that all patients identified as having moderate to end-stage (stages 3–5) kidney disease in whom a GBMCA is to be administered provide informed consent when practical, which includes a review of known risks and benefits as well as the possible availability of alternative imaging methods, if any.

As noted above, early data suggest that elevated levels of phosphate, iron, zinc, or copper might serve as efficient competitors for the "attention" of the chelate molecule [46]. These might therefore result

in increased levels of free gadolinium (Gd³⁺) ion in the patient, which might in turn increase the potential of the patient to develop NSF. Other cations such as lanthanum, now used as lanthanum carbonate (Fosrenol) for phosphorus binding in end-stage renal disease patients, could also present similar transmetallation and free gadolinium concerns. A history of multiple prior GBMCA administrations or hepatorenal disease also seems to be associated with an increased incidence of subsequent development of NSF. The existence of acidosis or active inflammatory and/or thrombotic processes as possible increased risk factors has been entertained but has not been reproducibly established at this point. This information should be taken into account during the risk-benefit assessment of each individual patient.

For administration of GBMCAs to patients on hemodialysis, the patient is to be transported to hemodialysis immediately upon the termination of the MR imaging examination. Arrangements should be made with the treating dialysis centers to provide them with as much notice as possible prior to the arrival of that patient for hemodialysis. It is recommended that hemodialysis be initiated no later than 2 hours following the administration of the GBMCA. This applies equally to emergent or urgent gadolinium chelate administration to these patients and to inpatients as well as outpatients. An additional hemodialysis session should be considered within 24 hours of this first contrast-enhanced treatment session for the reasons noted above.

For administration to patients on continuous ambulatory peritoneal dialysis (CAPD) or continuous cycling peritoneal dialysis (CCPD) (also known as automated peritoneal dialysis, or APD), there appears to be strong reason to hesitate to administer these agents. As noted above, this process of dialysis seems to be relatively ineffective at clearing the gadolinium from the body. Thus, special caution should be exercised when deciding whether a peritoneal dialysis patient should receive gadolinium-based MR contrast agents. If it is decided that they should be administered such agents, administration of the lowest reasonable dose is strongly recommended. In the past, it had been recommended that the patient avoid periods of a dry abdomen (i.e., no dialysate in the peritoneal cavity) and that the patient be advised to begin additional dialysis self-treatments or CCPD treatments immediately upon the termination of the MR examination in which the GBMCA was administered. These suggestions seemed prudent, although the efficacy of these recommendations had not been established. However, in light of the near-total apparent ineffectiveness of peritoneal dialysis at clearing the gadolinium from the body, it may well be worth considering immediate initiation of hemodialysis in peritoneal dialysis patients who receive even a low dose of a GBMCA, or not administering the agent if clinically feasible. Investigations are ongoing at this time to attempt to assess prevalence rates of NSF in peritoneal dialysis versus hemodialysis patients, although at this time it is too early for definitive conclusions.

Historically, as a result of the high atomic number associated with GBMCAs, these agents have occasionally been administered to (especially renal failure) patients in an off-label manner for such X-ray-based diagnostic tests as conventional angiography (including access angiography and fistulography) and even CT scanning. The rationale behind this practice was to avoid the administration of iodinated contrast agents to these patients and to decrease the incidence or likelihood of the development of contrast-induced nephropathy. In an attempt to prevent inadvertent GBMCA administration to renal disease patients by nonradiologists (who may at this point still not be fully aware of the issues and risks associated with GBMCAs), for now it is thought pr-

dent to ensure that all GBMCAs are to be administered only by radiologists. If there is a request for a GBMCA to be administered by a nonradiologist to a patient for an off-label use, such as intraarterial administration for vascular assessment in renal failure patients, this must be made in the form of a written order. All such requests must be prospectively reviewed and approved by either a radiologist or a pharmacist knowledgeable in the issues raised above, a risk-benefit assessment should be prospectively performed, and, where practical, informed consent should be provided by the patient.

For patients in whom a diagnosis of NSF has already been established, it might be appropriate to consider avoiding entirely any administration of a gadolinium-based MR contrast agent.

For patients not already on hemodialysis, the FDA's December 22, 2006 advisory [47] notwithstanding, the decision to initiate hemodialysis following gadolinium administration should not be taken lightly. The vast majority of NSF cases developed in patients with severe or end-stage renal disease, and most were already dialysis patients. The numbers of patients with moderate, as opposed to severe or end-stage, renal disease who have been diagnosed with NSF is exceedingly small, if they exist at all. At this time, it seems reasonable to assume that as the renal function/GFR decreases from 60 mL/min/1.73 m² through 30 mL/min/1.73 m², 15 mL/min/1.73 m², and below, the greater the concern and the greater the likelihood of subsequent NSF development. Therefore, we think that at the present time insufficient data exist to advise consideration for hemodialysis in this population of patients with moderate chronic kidney disease (stage 3) in the same manner or same perceived risk as those with severe or end-stage renal disease (stages 4 and 5). The risks of initiating hemodialysis must be seriously weighed against those of developing NSF in each particular case before a decision is made one way or the other. Finally, withholding clinically indicated GBMCAs can also be associated with its own risks, which should be considered in the decision-making process for all patients with kidney disease.

Should a new diagnosis of NSF be made, it is recommended that the FDA be notified through their MedWatch program (<http://www.fda.gov/medwatch/>) [11] or by phone (1-800-FDA-1088), and that the international NSF registry at Yale University be notified as well (<http://www.icnfd.org>) [39] to ensure that each database is kept as current as possible on this rapidly changing environment.

In the weeks and months to come, it is anticipated that there will be much further study of this issue, and that more information will be forthcoming that will hopefully shed more light on this important issue [56].

M. Patients in Whom There Are or May Be Intracranial Aneurysm Clips

1. In the event that it is unclear whether a patient does or does not have an aneurysm clip in place, plain films should be obtained. Alternatively, if available, any cranial plain films, CT, or MR examination that may have taken place in the recent past (i.e., subsequent to the suspected surgical date) should be reviewed to assess for a possible intracranial aneurysm clip.
2. In the event that a patient is identified to have an intracranial aneurysm clip in place, the MR examination should not be performed until it can be documented that the type of aneurysm clip within that patient is MR safe or MR conditional. All documentation of types of implanted clips, dates, etc., *must* be in writing and signed by a licensed physician. Phone or verbal histories and histories pro-

vided by a nonphysician are not acceptable. Fax copies of operative reports, physician statements, etc. are acceptable as long as a legible physician signature accompanies the requisite documentation. A written history of the clip itself having been appropriately tested for ferromagnetic properties (and description of the testing methodology used) prior to implantation by the operating surgeon is also considered acceptable if the testing follows the deflection test methodology established by ASTM International.

3. All implanted intracranial aneurysm clips that are documented in writing to be composed of titanium (either the commercially pure or the titanium alloy types) can be accepted for scanning without any other testing.
4. All nontitanium intracranial aneurysm clips manufactured in 1995 or later for which the manufacturer's product labeling continues to claim MR compatibility may be accepted for MR scanning without further testing.
5. Clips manufactured prior to 1995 require either pretesting (according to the ASTM deflection test methodology) prior to implantation or individual review of previous MRI of the clip or brain in that particular case, if available. By assessing the size of the artifact associated with the clip relative to the static field strength on which it was studied, the sequence type, and the MRI parameters selected, an opinion may be issued by one of the site's level 2 MR attending radiologists as to whether the clip demonstrates significant ferromagnetic properties or not. Access to the MR scanner would then be based on that opinion.
6. A patient with an aneurysm clip (or other implant) may have safely undergone a prior MR examination at any given static magnetic field strength. This fact in and of itself is not sufficient evidence of the implant's MR safety and should not solely be relied upon to determine the MR safety or compatibility status of that aneurysm clip (or other implant).

Variations in static magnetic field strength, static magnetic field spatial gradient, orientation of the aneurysm clip (or other implant) to the static magnetic field or static field gradient, rate of motion through the spatial static field gradient, etc., are all variables that are virtually impossible to control or reproduce. These variables may not have resulted in an adverse event in one circumstance but may result in significant injury or death on a subsequent exposure. For example, a patient who went blind from interactions between the metallic foreign body in his retina and the spatial static fields of the MR scanner entered the magnet and underwent the entire MR examination without difficulty; he went blind only on exiting the MR scanner at the completion of the examination.

7. Barring availability of either pretesting or prior MRI data of the clip in question, a risk-benefit assessment and review must be performed in each case individually. Further, for patients with intracranial clips with no available ferromagnetic or imaging data, should the risk-benefit ratio favor the performance of the MR study, the patient or guardian should provide written informed consent that includes death as a potential risk of the MRI procedure before that patient is permitted to undergo an MR examination.

N. Patients in Whom There Are or May Be Cardiac Pacemakers or Implantable Cardioverter Defibrillators

It is recommended that the presence of implanted cardiac pacemakers or implantable cardioverter defibrillators (ICDs) be considered a rel-

ative contraindication for MRI. MRI of patients with pacemakers and ICDs ("device patients") is *not* routine. Should an MRI be considered, it should be done on a case-by-case and site-by-site basis, and only if the site is staffed with individuals with the appropriate radiology and cardiology knowledge and expertise on hand. As of this writing, no cardiac pacing and/or defibrillating devices are labeled safe or conditionally safe for MRI scanning. Pacemaker and/or ICD leads may also present a hazard in the absence of any implant connected to them.

The protective circuitry of pacemakers and ICDs and their resistance to electromagnetic interference (EMI) has steadily improved over the years. Therefore, recently marketed ("modern") devices may be safer in the MRI environment. However, the committee eschews the term "modern" when referring to a particular device, recognizing that all devices currently marketed contain legacy components that may or may not be resistant to the forces and EMI present in the MRI suite. Future devices, unless appropriately tested and labeled as such, should not be regarded as safe for MRI simply because they are "modern" or recently manufactured.

Unexpected programming changes, inhibition of pacemaker output, failure to pace, transient asynchronous pacing, rapid cardiac pacing, the induction of ventricular fibrillation, heating of the tissue adjacent to the pacing or ICD system, early battery depletion, and outright device failure requiring replacement may all occur during MRI of patients with pacemakers or ICDs. The committee notes that multiple deaths have occurred under poorly and incompletely characterized circumstances when device patients underwent MRI. These deaths may have occurred as a result of pacemaker inhibition, failure to capture or device failure (resulting in prolonged asystole), and/or rapid cardiac pacing or asynchronous pacing (resulting in the initiation of ventricular tachycardia or fibrillation).

Ideally, the nonemergent patient should be apprised of the risks associated with the procedure and should provide prospective written informed consent prior to the initiation of MRI. While the majority of reported deliberate scans of device patients have proceeded without mishap when appropriate precautions were taken, there may be underreporting of adverse events, including deaths [58]. Thus, assignment of a risk-benefit ratio to the performance of MRI in a device patient is difficult. While the risk may be low, device patients who are considered for MRI should be advised that life-threatening arrhythmias might occur during MRI and serious device malfunction might occur, requiring replacement of the device.

Should any MRI examination be contemplated for a patient with an implanted pacemaker or ICD, it is recommended that radiology and cardiology personnel and a fully stocked crash cart be readily available throughout the procedure in case a significant arrhythmia develops during the examination that does not terminate with the cessation of the MR study. The cardiologist should be familiar with the patient's arrhythmia history and the implanted device. A programmer that can be used to adjust the device as necessary should be readily available. All such patients should be actively monitored for cardiac and respiratory function throughout the examination. At a minimum, ECG and pulse oximetry should be used. At the conclusion of the examination, the cardiologist should examine the device to confirm that the function is consistent with its preexamination state. Follow-up should include a check of the patient's device at a time remote (1–6 weeks) after the scan to confirm appropriate function.

Should an MRI (or entry into the magnet area) be performed inadvertently on a patient with a pacemaker or ICD, the patient's car-

diologist should be contacted before the patient's discharge from the MRI suite. The importance of examination of the device prior to the patient's leaving the MRI suite cannot be overstated.

O. Site Emergency Preparedness

There are many factors to consider when attempting to ensure that an MR imaging facility is adequately and appropriately prepared to handle any of several types of emergencies that might impact MR scanners of varied design types. Appendix 4 addresses these issues in some detail and provides specific guidelines to help anticipate and safeguard sites from some of the more common emergencies and disasters that might affect MR imaging facilities.

References

1. Chaljub G, Kramer LA, Johnson RF III, Singh H, Crow WN. Projectile cylinder accidents resulting from the presence of ferromagnetic nitrous oxide or oxygen tanks in the MR suite. *AJR* 2001; 177:27–30
2. ECRI hazard report. Patient death illustrates the importance of adhering to safety precautions in magnetic resonance environments. *Health Devices* 2001; 30:311–314
3. Kanal E, Borgstede JP, Barkovich AJ, et al. American College of Radiology white paper on MR safety. *AJR* 2002; 178:1335–1347
4. Kanal E, Borgstede JP, Barkovich AJ, et al. American College of Radiology white paper on MR safety: 2004 update and revisions. *AJR* 2004; 182:1111–1114
5. Shellock FG, Kanal E. Policies, guidelines and recommendations for MR imaging safety and patient management. The SMRI Safety Committee. *J Magn Reson Imaging* 1991; 1:97–101
6. Kanal E, Shellock FG. Policies, guidelines and recommendations for MR imaging safety and patient management. The SMRI Safety Committee. *J Magn Reson Imaging* 1992; 2:247–248
7. Shellock FG, Kanal E. Guidelines and recommendations for MR imaging safety and patient management. III. Questionnaire for screening patients before MR procedures. The SMRI Safety Committee. *J Magn Reson Imaging* 1994; 4:749–751
8. American College of Radiology. ACR standard for magnetic resonance imaging monograph: safety and sedation. In: *Practice guidelines and technical standards* 1996. Reston, VA: American College of Radiology, 1996
9. American College of Radiology. ACR standard for performing and interpreting magnetic resonance imaging (MRI). In: *Practice guidelines and technical standards* 2000. Reston, VA: American College of Radiology, 2000:429–433
10. Medical Device Reporting (MDR) page. Center for Devices and Radiological Health, Food and Drug Administration Website. Available at: www.fda.gov/cdrh/devadvice/351.html. Accessed February 6, 2007
11. U.S. Food and Drug Administration. MedWatch: The FDA Safety Information and Adverse Event Reporting Program. Available at: www.fda.gov/medwatch/. Accessed February 12, 2007
12. Jarvik JG, Ramsey S. Radiographic screening for orbital foreign bodies prior to MR imaging: is it worth it? *Am J Neuroradiol* 2000; 21:245–247
13. Seidenwurm DJ, McDonnell CH III, Raghaven N, Breslau J. Cost utility analysis of radiographic screening for an orbital foreign body before MR imaging. *Am J Neuroradiol* 2000; 21:426–433
14. Kanal E, Gillen J, Evans JA, Savitz DA, Shellock FG. Survey of reproductive health among female MR workers. *Radiology* 1993; 187:395–399
15. American Academy of Pediatrics Committee on Drugs. Guidelines for monitoring and management of pediatric patients during and after sedation for diagnostic and therapeutic procedures. *Pediatrics* 1992; 89(6 Pt 1):1110–1115
16. American Academy of Pediatrics Committee on Drugs. Guidelines for monitoring

- and management of pediatric patients during and after sedation for diagnostic and therapeutic procedures: addendum. *Pediatrics* 2002; 110:836–838
17. American Society of Anesthesiologists. *Updated practice guidelines for sedation and analgesia by non-anesthesiologists*. Park Ridge, IL: American Society of Anesthesiologists, 2001
 18. Joint Commission on the Accreditation of Healthcare Organizations. *Standards and intents for sedation and anesthesia care: comprehensive accreditation manual for hospitals*. Report no. TX 2-2. Chicago, IL: JCAHO, 2002
 19. Bryan YF, Templeton TW, Nick TG, Szafran M, Tung A. Brain magnetic resonance imaging increases core body temperature in sedated children. *Anesth Analg* 2006; 102:1674–1679
 20. Gandhi OP, Chen XB. Specific absorption rates and induced current densities for an anatomy-based model of the human for exposure to time varying magnetic fields of MRI. *Magn Reson Med* 1999; 41:816–823
 21. Konings MK, Bartels LW, Smits HFM, et al. Heating around intravascular guidewires by resonating RF waves. *J Magn Reson Imaging* 2000; 12:79–85
 22. Shellock FG. *Pocket guide to MR procedures and metallic objects: update 2001*. Philadelphia, PA: Lippincott Williams & Wilkins, 2001:147–148
 23. International Electrotechnical Commission. *Medical electrical equipment. Part 2. Particular requirements for the safety of magnetic resonance equipment for medical diagnosis. IEC 60601-2-33, ed. 2.0*. Geneva, Switzerland: IEC, 2002. (In the United States, copies of this standard can be obtained from the American National Standards Institute (ANSI), 11 W 42nd St., New York, NY 10036.)
 24. American College of Radiology. ACR practice guideline for adult sedation/analgesia. In: *Practice guidelines and technical standards 2000*. Reston, VA: American College of Radiology, 2000
 25. American College of Radiology. ACR practice guideline for pediatric sedation/analgesia. In: *Practice guidelines and technical standards 1998*. Reston, VA: American College of Radiology, 1998
 26. American Society of Anesthesiologists. *Guidelines for non-operating room anesthetizing locations*. Park Ridge, IL: American Society of Anesthesiologists, 2001
 27. American Society of Anesthesiologists. *Standards for basic anesthetic monitoring*. Park Ridge, IL: American Society of Anesthesiologists, 1986, 1998
 28. American Society of Anesthesiologists. *Standards for post anesthesia care*. Park Ridge, IL: American Society of Anesthesiologists, 1988, 1994
 29. Joint Commission on Accreditation of Healthcare Organizations. *Standards and intents for sedation and anesthesia care: comprehensive accreditation manual for hospitals*. Report TX. 2-2.4.1. Chicago, IL: JCAHO, 2001
 30. American College of Radiology. *ACR digest of council actions*. Reston, VA: ACR, 1999:126 (Res. 1-H, 1987, 1997)
 31. American College of Radiology Committee on Drugs and Contrast Media. *Manual on contrast media*, 5.0 ed. Reston, VA: American College of Radiology, 2004
 32. Grobner T. Gadolinium: a specific trigger for the development of nephrogenic fibrosing dermopathy and nephrogenic systemic fibrosis? *Nephrol Dial Transplant* 2006; 21:1104–1108; erratum in *Nephrol Dial Transplant* 2006; 21:1745
 33. Investigation of the safety of MRI contrast medium Omniscan, May 29, 2006. Danish Medicines Agency Website. Available at: www.dkma.dk/1024/visUKLSArtikel.asp?artikelID=8931. Accessed January 25, 2007; last updated January 24, 2007
 34. Cowper SE, Robin HS, Steinberg SM, Su LD, Gupta S, LeBoit PE. Scleromyxedema-like cutaneous diseases in renal-dialysis patients. *Lancet* 2000; 356:1000–1001
 35. Jain SM, Wesson S, Hassanein A, et al. Nephrogenic fibrosing dermopathy in pediatric patients. *Pediatr Nephrol* 2004; 19:467–470
 36. Dharmidharka VR, Wesson SK, Fennell RS. Gadolinium and nephrogenic fibrosing dermopathy in pediatric patients. (letter) *Pediatr Nephrol* 2006 Dec 16; [Epub ahead of print]
 37. Broome DR, Girguis MS, Baron PW, Cottrell AC, Kjellin I, Kirk GA. Gadodiamide-associated nephrogenic systemic fibrosis: why radiologists should be concerned. *AJR* 2007; 188:586–592
 38. Maloo M, Abt P, Younan D, et al. Nephrogenic systemic fibrosis among liver transplant recipients: a single institution experience and topic update. *Am J Transplant* 2006; 6:2212–2217
 39. Cowper SE. Nephrogenic fibrosing dermopathy [NFD/NSF Website], 2001–2007. Available at: www.icnldr.org. Accessed January 24, 2007; site last updated January 24, 2007
 40. Marckmann P, Skov L, Rossen K, et al. Nephrogenic systemic fibrosis: suspected causative role of gadodiamide used for contrast-enhanced magnetic resonance imaging. *J Am Soc Nephrol* 2006; 17:2359–2362
 41. Thomsen HS. Nephrogenic systemic fibrosis: a serious later adverse reaction to gadodiamide. *Eur Radiol* 2006; 16:2619–2621
 42. Khurana A, Runge VM, Narayanan M, Greene JF Jr, Nickel AE. Nephrogenic systemic fibrosis: a review of 6 cases temporally related to gadodiamide injection (Omniscan). *Invest Radiol* 2007; 42:139–145
 43. White GW, Gibby WA, Tweedle MF. Comparison of Gd (DTPA-BMA) (Omniscan) versus Gd (HP-DO3A) (ProHance) relative to gadolinium retention in human bone tissue by inductively coupled plasma mass spectroscopy. *Invest Radiol* 2006; 41:272–278
 44. High WQ, Ayers RA, Chandler J, Zito G, Cowper SE. Gadolinium is detectable within the tissues of patients with nephrogenic systemic fibrosis. *J Am Acad Dermatol* 2007; 56:21–26
 45. Boyd AS, Zic JA, Abraham JL. Gadolinium deposition in nephrogenic fibrosing dermopathy. *J Am Acad Dermatol* 2007; 56:27–30
 46. Idee J-M, Port M, Raynal I, Schaefer M, Greneur SL, Corot C. Clinical and biological consequences of transmetallation induced by contrast agents for magnetic resonance imaging: a review. *Fundam Clin Pharmacol* 2006; 20:563–576
 47. Food and Drug Administration Center for Drug Evaluation and Research. Public health advisory. Update on magnetic resonance imaging (MRI) contrast agents containing gadolinium and nephrogenic fibrosing dermopathy. Available at: www.fda.gov/cder/drug/advisory/gadolinium_agents_20061222.htm. Accessed February 20, 2007; last updated December 22, 2006
 48. U.S. Food and Drug Administration Website. Public health advisory. Gadolinium-containing contrast agents for magnetic resonance imaging (MRI): Omniscan, OptiMARK, Magnevist, ProHance, and MultiHance. Available at: http://www.fda.gov/cder/drug/advisory/gadolinium_agents.htm. Accessed August 10, 2006
 49. Okada S, Katagiri K, Kumazaki T, Yokoyama H. Safety of gadolinium contrast agent in hemodialysis patients. *Acta Radiol* 2001; 42:339–341
 50. Joffe P, Thomsen HS, Meusel M. Pharmacokinetics of gadodiamide injection in patients with severe renal insufficiency and patients undergoing hemodialysis or continuous ambulatory peritoneal dialysis. *Acad Radiol* 1998; 5:491–502
 51. Ueda J, Furukawa T, Higashino K, et al. Permeability of iodinated and MR contrast media through two types of hemodialysis membrane. *Eur J Radiol* 1999; 31:76–80
 52. The National Kidney Foundation Kidney disease outcomes quality initiative. National Kidney Foundation Website. Available at: <http://www.kidney.org/Professionals/kdoqi/>. Accessed February 12, 2007
 53. Centers for Disease Control and Prevention, National Center for Health Statistics Website. Third National Health and Nutrition Examination Survey (NHANES III), 1988–94. Available at: www.cdc.gov/nchs/products/elec_prods/subject/nhanes3.htm. Accessed February 8, 2007
 54. [No authors indicated] Part 4. Definition and classification of stages of chronic kidney disease. *Am J Kidney Dis* 2002; 39[suppl 1]:S46–S75; doi:10.1053/ajkd.2002.30943
 55. Kuo PH, Kanal E, Abu-Alfa AK, Cowper SE. Gadolinium-based MR contrast agents and nephrogenic systemic fibrosis. *Radiology* 2007 Jan 9 [Epub ahead of print]. Available at: <http://radiology.rsnajnl.org/cgi/content/>

short/2423061640v1?rss=1. Accessed February 20, 2007

56. Medicines and Healthcare products Regulatory Agency Website. Nephrogenic systemic fibrosis (NSF) and gadolinium-containing MRI contrast agents: 7 February 2007. Available at: http://www.mhra.gov.uk/home/idcplg?IdcService=SS_GET_PAGE&useSecondary=true&ssDocName=CON2030229&ssTargetNodeId=221. Accessed February 20, 2007
57. Thomsen HS. European Society of Urogenital Radiology guidelines on contrast media application. *Curr Opin Urol* 2007; 17:70–76
58. Bartsch C, Irnich W, Risse M, Weiler G. Unexpected sudden death of pacemaker patients during or shortly after magnetic resonance imaging (MRI). In: *XIX Congress of International Academy of Legal Medicine*. IALM: Milan, Italy, 3–6 September 2003:174 (Abstract No. 114)
1. Luechinger R, Zeijlemaker VA, Pederson EM, et al. In vivo heating of pacemaker leads during magnetic resonance imaging. *Eur Heart J* 2005; 26:376–383
2. Nitz WR, Brinker G, Diehl D, Frese G. Specific absorption rate as a poor indicator of magnetic resonance-related implant heating. *Invest Radiol* 2005; 40:773–776
3. Nguyen UD, Brown SJ, Chang IA, et al. Numerical evaluation of heating of the human head due to magnetic resonance imaging. *IEEE Transactions on Biomedical Engineering* 2004; 51:1301–1309
4. Nyenhuis JA, Park SM, Kamondetdacha R, et al. MRI and implanted medical devices: basic interactions with an emphasis on heating. *IEEE Transactions on Device and Materials Reliability* 2005; 5:467–480
5. Shellock FG, Fieno DS, Thomson LJ, Talavage TM, Berman DS. Cardiac pacemaker: in vitro assessment at 1.5 T. *Am Heart J* 2006; 2:436–443
6. Yeung CJ, Susil RC, Atalar E. RF safety of wires in interventional MRI: using a safety index. *Magn Reson Med* 2002; 47:187–193
7. Shellock FG, Begnaud J, Inman MD. Vagus nerve stimulation therapy system: in vitro evaluation of magnetic resonance imaging-related heating and function at 1.5 and 3 Tesla. *Intermodulation* 2006; 9:204–213
8. Irnich W, Irnich B, Bartsch C, et al. Do we need pacemakers resistant to magnetic resonance imaging? *Europace* 2005; 7:353–365
9. Loewy J, Loewy A, Kendall EJ. Reconsideration of pacemakers and MR imaging. *RadioGraphics* 2004; 24:1257–1268
10. Luechinger R, Duru F, Candinas R, Boesiger P. Safety considerations for magnetic resonance imaging of pacemaker and ICD patients [in German]. *Herzschrittmachertherapie und Elektrophysiologie* 2004; 15:73–81
11. ASTM International. *Terminology of symbols and definitions relating to magnetic testing*, A 340. West Conshohocken, PA: ASTM International, 2003:1–16
12. ASTM International. *Specification for the requirements and disclosure of self-closing aneurysm clips*, F 1542. West Conshohocken, PA: ASTM International, 2000:1–3
13. ASTM International. *Test method for measurement of magnetically induced displacement force on medical devices in the magnetic resonance environment*, F 2052. West Conshohocken, PA: ASTM International, 2002: 1–6
14. ASTM International. *Test method for evaluation of MR image artifacts from passive implants*, F 2119. West Conshohocken, PA: ASTM International, 2001:1–3
15. ASTM International. *Test method for measurement of radio frequency induced heating near passive implants during magnetic resonance imaging*, F 2182. West Conshohocken, PA: ASTM International, 2002:1–8
16. ASTM International. *Standard test method for measurement of magnetically induced torque on passive implants in the magnetic resonance environment*, F 2213. West Conshohocken, PA: ASTM International, 2004:1–7

Additional Sources

Appendices 1–4 appear on the following pages

FOR YOUR INFORMATION

See the accompanying CME/SAM article in this month's issue of *AJR Integrative Imaging* included with your *AJR*.

FOR YOUR INFORMATION

The reader's attention is directed to the commentary on this article, which appears on page 1446.

APPENDIX 1: Personnel and Zone Definitions

PERSONNEL DEFINITIONS

Non-MR Personnel

Patients, visitors, or facility staff who do not meet the criteria of level 1 or level 2 MR personnel will be referred to as non-MR personnel. Specifically, non-MR personnel will be the terminology used to refer to any individual or group who has not within the previous 12 months undergone the designated formal training in MR safety issues defined by the MR safety director of that installation.

Level 1 MR Personnel

Individuals who have passed minimal safety educational efforts to ensure their own safety as they work in Zone III will be referred to as level 1 MR personnel (e.g., MRI department office staff and patient aides).

Level 2 MR Personnel

Individuals who have been more extensively trained and educated in the broader aspects of MR safety issues, including issues related to the potential for thermal loading or burns and direct neuromuscular excitation from rapidly changing gradients, will be referred to as level 2 MR personnel (e.g., MRI technologists, radiologists, and radiology department nursing staff).

ZONE DEFINITIONS

Zone I

This region includes all areas that are freely accessible to the general public. This area is typically outside the MR environment itself and is the area through which patients, health care personnel, and other employees of the MR site access the MR environment.

Zone II

This area is the interface between the publicly accessible uncontrolled Zone I and the strictly controlled Zone III (see below). Typically, the patients are greeted in Zone II and are not free to move throughout Zone II at will, but rather are under the supervision of MR personnel. It is in Zone II that the answers to MR screening questions, patient histories, medical insurance questions, etc. are typically obtained.

Zone III

This area is the region in which free access by unscreened non-MR personnel or ferromagnetic objects or equipment can result in serious injury or death as a result of interactions between the individuals or equipment and the particular environment of the MR scanner. These interactions include, but are not limited to, those with the MR scanner's static and time-varying magnetic fields. All access to Zone III is to be strictly restricted, with access to regions within it (including Zone IV, see below) controlled by, and entirely under the supervision of, MR personnel.

Zone IV

This area is synonymous with the MR scanner magnet room itself. Zone IV, by definition, will always be located within Zone III as it is the MR magnet and its associated magnetic field which generates the existence of Zone III.

Non-MR personnel should be accompanied by, or under the immediate supervision of and visual contact with, one specifically identified level 2 MR person for the entirety of their duration within Zone III or Zone IV restricted regions.

Levels 1 and 2 MR personnel may move freely about all zones.

APPENDIX 2: MR Facility Safety Design Guidelines

The goal of MR safety is to prevent harm to patients, though an MR facility cannot simply adopt one or two interventions and hope to successfully attain this objective. According to safety and human factors engineering principles, multiple safety strategies must be adopted to be effective. This approach is sometimes termed "defense in depth." The safety strategies outlined in the main body of this Guidance Document for Safe MR Practices include, for instance, policies that restrict personnel access, specialized training and drills for MR personnel, and warning labels for devices to be brought into Zone IV regions.

Along with these people-oriented strategies of policies and training, organizations need also to adopt the strategies of safety-oriented architectural and interior design. These design elements can support the other safety strategies by making them easier or more obvious to follow. The architectural enhancements described herein add one or more strong barriers to enhance "defense in depth."

This appendix includes descriptions of architectural and interior design recommendations organized around the many MR suite functional areas. Note that a facility's design can encourage safety and best practices by improving the flow of patients, various health care per-

sonnel, and equipment and devices, and not just by preventing MR unsafe items from becoming missiles, or screening out patients with hazardous implanted devices.

Placing design elements strategically in a suite layout such that the element supports best-practice work flow patterns will increase compliance with safer practices. For example, having a private area for patient screening interviews will make it more likely the patients will disclose sensitive types of implants. Another example of designing for safety is to include dedicated space and temporary storage for MR Unsafe equipment (e.g., ferromagnetic IV poles, transport stretchers) out of direct sight and away from people flow patterns.

Effective and safe MRI suites must balance the technical demands of the MR equipment with local and state building codes, standards of accrediting bodies, clinical and patient population needs, payor requirements, and a collage of civil requirements from the Health Insurance Portability and Accountability Act (HIPAA) to the Americans with Disabilities Act (ADA).

In an effort to better match appropriate facility design guidelines with levels of patient acuity and care, the ACR MR Safety Committee

is currently developing level designations for MRI facilities in conjunction with the efforts of committees from other societies and organizations. These will address customization of requirements for sites with varying anticipated patient care sedation, anesthesia, and/or interventional activities.

While it would be desirable to provide a universal MRI suite safety design, the variables are too numerous to adequately address in a single template. The following MRI Facility Safety Design Guidelines provide information in support of planning, design, and construction of MR facilities, including updates to existing MR facilities, which enhance the safety of patients, visitors, and staff. This information is intended to supplement and expand upon patient safety guidance provided throughout the ACR Guidance Document for MR Safe Practices.

1. MR Equipment Vendor Templates

Design templates provided by MR equipment manufacturers are invaluable in developing suites that meet the minimum technical siting requirements for the specific equipment. Vendor design templates, however, typically depict only the control and equipment rooms, in addition to the magnet room, Zone IV.

Patient/family waiting, interview areas, physical screening/changing areas, access controls, storage, crash carts, induction, medical gas services, holding areas for patients after screening, infection control provisions, and interventional applications, among many other issues, are not addressed in typical vendor-provided drawings. These issues are left to facility owners, operators, and their design professionals to resolve. The guidance which follows is designed to address many of these issues which directly impact safety within the MR suite.

2. Patient Interview/Clinical Screening Areas (Zone II)

Reviewing the patient Safety Screening Form and the MR Hazard Checklist requires discussing confidential personal information. To facilitate full and complete patient disclosure of their medical history, this clinical screening should be conducted in an area which provides auditory and visual privacy for the patient. Facilities should prospectively plan for electronic patient medical records, which are useful in clinical screening, and should provide access to records in the MR suite in support of clinical patient screening.

Clinical screening of inpatients may be completed in the patient room for hospital-based MR facilities. However, all screenings are to be double-checked and verified by appropriately trained MR personnel before MR examination.

3. Physical Screening and Patient Changing/Gowning Rooms (Zone II)

All persons and objects entering Zone III should be physically screened for the presence of ferromagnetic materials which, irrespective of size, can become threats in proximity to the MR scanner. A location should be provided for patients in which they may change out of their street clothes and into a facility-provided gown or scrubs, if or as deemed appropriate. For those facilities that either do not provide space for, or do not require, patient changing, the facility must provide alternative means of identifying and removing items that the patient may have brought with them that might pose threats in the MR environment.

A high-strength handheld magnet is a recommended tool to evaluate the gross magnetic characteristics of objects of unknown composition. Magnetic strength for these permanent magnets falls off quickly as one moves away from the face of the magnet. Thus, these may not demon-

strate attraction for ferromagnetic components which are not superficially located or cannot for whatever reason be brought into close proximity with the surface of this handheld magnet.

Ferromagnetic detection systems have been demonstrated to be highly effective as a quality assurance tool, verifying the successful screening and identifying ferromagnetic objects which were not discovered by conventional screening methods. It is recommended that new facility construction anticipate the use of ferromagnetic detection screening in Zone II and provide for installation of the devices in a location which facilitates use and throughput. Many current ferromagnetic detection devices are capable of being positioned within Zone III, even at the door to the magnet room; however, the recommended use of ferromagnetic detection is to verify the screening of patients before they pass through the controlled point of access into Zone III.

Physical screening of patients should consist of removal of all jewelry, metallic or ferromagnetic objects, onplants, and prostheses (as indicated by manufacturer's conditional use requirements and physician instructions) and either having patients change out of their street clothes into facility-provided gowns or scrubs or thorough screening of street clothes, including identifying the contents of pockets and the composition of metallic fibers, fasteners, and reinforcing.

4. Transfer Area and Ferrous Quarantine Storage (Zone II)

Patients arriving with wheelchairs, walkers, portable oxygen, and other appliances that may be unsafe in the MR environment should be provided by the facility with appropriate MR safe or MR conditional appliances. An area should be provided to transfer the patient from unsafe appliances to ones appropriate to the MR environment. Unsafe appliances brought by the patient should be secured in a "ferrous quarantine" storage area, distinct from storage areas for MR safe and MR conditional equipment, and ideally locked out of sight. Patient belongings should be retrieved from the ferrous quarantine only when discharging the patient to whom the objects belong.

5. Access Control (Zone III and Zone IV)

Means of physically securing and restricting access to Zone III from all adjacent areas must be provided. Independent access into Zone III must be limited to only appropriately trained MR personnel.

6. Patient Holding (Zone III)

Depending upon facility capacity and patient volume, it may be advisable to provide a postscreening patient holding area. Zone III holding areas should be equipped and appointed to prevent patient exit and subsequent reentry. This will help prevent the inadvertent—or even intentional—introduction of unscreened objects and personnel.

Many multitechnique radiology facilities combine patient holding and/or induction areas for patients undergoing different types of imaging examinations. This presents safety challenges when, for example, patients scheduled to undergo CT are held in a patient holding area shared by postscreening MR patients. As CT patients would not typically be screened for MR contraindications or ferrous materials, this poses risks to both the CT patient with a contraindicated implant and to those in the MRI Zone IV should an unscreened individual inadvertently enter with a ferrous object or implant.

Unless all persons in patient holding areas used for postscreening MR patients are screened for MRI, the practice of shared patient holding areas between MR and other techniques is discouraged. Ultimately it is the re-

sponsibility of trained MR staff to verify the screening of any commingled patient prior to permitting them to enter Zone III and Zone IV.

In all MR facilities, Zone III is required to be secured and access limited to only MR personnel and successfully MR prescreened non-MR personnel accompanied by MR personnel. Ideally, facilities should be designed so that patients undergoing other techniques are not commingled with postscreening MR patients.

7. Lines of Sight and Situational Awareness (Zone III)

Trained MR personnel are arguably the single greatest safety resource of MR facilities. These individuals should be afforded visual control over all persons entering or exiting Zone III or Zone IV. The technologist seated at the MR operator console should therefore be able to view not only the patient in the MR scanner but also the approach and entrance into Zone IV. When practical, suites should also be prospectively designed to provide a view from the MR operator's console to patient holding areas. If this cannot be satisfactorily achieved by direct line of sight, remote video viewing devices are an acceptable substitute toward accomplishing this objective.

The technologist at the console should also be provided with a view to induction and recovery areas within the MR suite, as applicable.

8. Emergency Resuscitation Equipment (Zone II or Zone III)

Because of risks associated with contrast agents, sedation, anesthesia, and even the frail health of patients undergoing MR examinations, it is advised that each facility have appropriate provisions for stabilization and resuscitation of patients.

It is recommended that crash carts and emergency resuscitation equipment be stored in a readily accessible area in either Zone II or Zone III. This emergency resuscitation equipment is to be appropriately labeled and also tested and verified as safe for usage in an MR environment for the anticipated conditions of usage.

MR facilities should maintain a supply of emergency medications to treat adverse reactions to administered contrast agents.

MR facilities providing care to patients who require clinical support during the MR examination should have emergency response equipment and personnel, trained in MR safety issues as well as trained to respond to anticipatable adverse events, readily available to respond to patient adverse events or distress in the MR arena.

9. Fringe Magnetic Field Hazards (Zone III)

For many MR system installations, magnetic fringe fields which project beyond the confines of the magnet room superimpose potential hazards on spaces which may be outside the MR suite, potentially on levels above or below the MR site and perhaps even outside the building. Facilities must identify all occupiable areas, including those outside the MR suite (including rooftops, storage areas, mechanical closets, etc.) which are exposed to potentially hazardous magnetic fringe field strengths. Areas of potential hazard must be clearly identified, and access to these areas must be restricted, just as they would be within the MR suite.

10. Cryogen Safety (Zone IV)

Liquid helium and liquid nitrogen represent the most commonly used cryogens in MR environments. The physical properties of these cryogenic liquids present significant potential safety hazards. If exposed to room air, these cryogenic liquids will rapidly boil off and expand into a gaseous state. This produces several potential safety concerns, including:

- Asphyxiation is a possibility as cryogenic gases replace oxygenated air.

- Frostbite may occur at the exceedingly low temperatures of these cryogenic liquids.
- Fire hazards can exist in the unlikely event of a quench, especially if some of the cryogenic gases escape into the magnet room/Zone IV.
- Hyperbaric pressure considerations within Zone IV can rarely exist in the unlikely event of a quench in which some of the cryogenic gases escape into the magnet room/Zone IV.

a. Cryogen Fills

Though contemporary superconducting magnets require cryogen refills only infrequently, there is still almost always the need to periodically bring hundreds of liters of liquid cryogen to the magnet. It is because of the risks to persons near the magnet and storage/transport dewars that transfill operations should be undertaken with great care and only by appropriately trained personnel.

- Dewars containing cryogenic liquids should never be stored inside an MRI facility or indeed any enclosed facility unless in a facility specifically designed to manage the associated pressure, temperature, and asphyxiation risks.
- A cryogen transfill should never be attempted by untrained personnel or even with any unnecessary personnel in attendance, including MR personnel staff and patients, within Zone IV.
- Cryogen transfills should only be performed with appropriate precautions in place to prevent pressure entrapment and asphyxiation.

b. Magnet Room Cryogen Safety

For most MRI systems, if the magnet quenches, the escaping cryogenic gases are ducted outside the building to an unoccupied discharge area. However, there have been documented failures of cryogen vent/quench pipe assemblies which have led to considerable quantities of cryogenic gases being inadvertently discharged into the magnet room/Zone IV. The thermal expansion of the cryogens, if released into the magnet room, can positively pressurize the magnet room and entrap persons inside until such time as the pressure is equalized.

The following recommended MRI suite design and construction elements reduce patient and staff risks in the unlikely event of a quench in which the cryogen vent pathway (quench pipe) ruptures or leaks into Zone IV:

- All magnet rooms/Zone IV regions for superconducting magnets should be provided with an emergency exhaust pathway. The emergency exhaust grille is to be located in the ceiling opposite the entrance to the magnet room (Zone IV) door. At this location, when activated in the unlikely event of a quench breach, the exhaust fan is positioned to draw the vaporous cloud of cryogenic gas away from the door providing exit from the magnet room.
- Many MR manufacturers are now requiring that magnet rooms for superconducting magnets also be provided with an additional form of passive pressure relief/pressure equalization to minimize the risks of positive-pressure entrapment. Designs for passive pressure relief mechanisms should follow design criteria similar to those of cryogen vent pathway and active exhaust, including discharge to a protected area, as described in section 10.c below.

Some MR facilities are constructed without open waveguides or glass observation windows to Zone IV regions. In these facilities, the potential risks of entrapment are even greater and may warrant an additional degree of attention in this regard.

While it can provide a degree of redundancy, it should be noted that, even with an exhaust fan, designing the door to Zone IV to swing out-

ward is not, by itself, an appropriate means of pressure relief. In a severe positive-pressure situation, unlatching an outward-swinging door might permit the door to burst open with tremendous pressure, potentially injuring person(s) opening the door. If employed as the only means of pressure equalization, an outward-swinging door may actually introduce new hazards to any staff person attempting to open the door to a pressurized magnet room from the outside.

Similarly, though it has proven effective in life-threatening situations, breaking a control window should not be advocated as a primary means of relieving/equalizing Zone IV pressure in a quench situation. It should be noted that the current construction of many RF-shielded observation windows is such that breaking the window would be very difficult, further diminishing that as a viable means of pressure relief.

Once provided with appropriate pressure equalization and emergency exhaust, magnet room door swing direction and design should be left to the discretion of a facility and their design professionals.

c. Cryogen Vent Pathway

Obstructions, inappropriate pipe materials, insufficient pipe caliber and/or length, or faulty connections in the length of the cryogen vent pathway can cause failure between the magnet and the point of discharge. An evaluation of the current cryogen vent piping/ducting assembly is recommended to help identify and correct potential weaknesses that could potentially fail in a quench. Facilities are advised to evaluate the design and inspect the construction of their cryogen vent system.

Because minimum design requirements for some cryogen vent systems have been revised by magnet system vendors, facilities should obtain current standards from the original equipment manufacturers to use in evaluating their cryogen vent assembly and not rely on original siting requirements.

Beyond the assessment of the current construction of the cryogen vent system, it is prudent for MRI facilities:

- To inspect cryogen vent systems at least annually, identifying stress or wear of pipe sections and couplings, loose fittings and supports, or signs of condensation or water within the cryogen vent pathway, which may indicate a blockage.
- Following any quench of a superconducting magnet, to conduct a thorough inspection of the cryogen vent system, including pipe sections, fittings, couplings, hangers, and clamps, prior to returning the magnet to service.

Because obstructions or occlusions of the cryogen vent can increase the likelihood of rupture in a quench event, facilities should ensure that:

- The discharge point has an appropriate weatherhead that prevents horizontal, wind-driven precipitation from entering, collecting, or freezing in the quench exhaust pipe.
- The discharge point is high enough off the roof or ground surface that snow or debris cannot enter or occlude the pipe.
- The discharge is covered by a material having sufficiently small openings to prevent birds or other animals from entering the quench pipe, while not occluding cryogenic gaseous egress in a quench situation.

Facilities that discover failings in any of these basic protections of the cryogen discharge point should immediately take additional steps to verify the patency of the cryogen vent and provide the minimum current discharge protections recommended by the original equipment manufacturer.

To protect persons from cryogen exposure at the point of discharge:

- At the point of cryogen discharge, a quench safety exclusion zone with a minimum clear radius of 25 ft (8 m) should be established and clearly marked with surface warnings and signage.
- The quench safety exclusion zone should be devoid of serviceable equipment, air intakes, operable windows, or unsecured doors that either require servicing or offer a pathway for cryogenic gasses to reenter the building.
- Persons who must enter this quench safety exclusion zone, including incidental maintenance personnel and contractors, should be permitted to do so only after receiving specific instruction on quench risks and response.

11. MR Conditional Devices (Zone IV)

The normal or safe operation of many medical devices designed for use in the MR environment may be disrupted by exposure to conditions exceeding the device's conditional rating threshold. It is advisable for MR facilities to identify the maximum conditional rating for static field and spatial gradient exposure for each MR Conditional device that may be brought into Zone IV. For prospective installations, it is recommended that the location of critical isogauss line(s) be identified for MR Conditional equipment and devices used within the MR suite and delineated on the floor and walls of the magnet room to aid in the positioning and safe and effective operation of said equipment.

All MR facilities should evaluate all MR Conditional patient monitoring, ventilators, medication pumps, anesthesia machines, monitoring devices, biopsy, and other devices and equipment which may be brought into the magnet room for magnetic field tolerances. Facilities should consider providing physical indications of critical gauss lines in the construction of the magnet room to promote the safe and effective use of MR Conditional equipment, as appropriate.

12. Infection Control (Zone IV)

Because of safety concerns regarding incidental personnel within the MR suite, restricting housekeeping and cleaning personnel from Zone III and/or Zone IV regions may give rise to concerns about the cleanliness of the MR suite. Magnet room finishes and construction details should be designed to facilitate cleaning by appropriately trained staff with nonmotorized equipment. Additionally, as the numbers of MR-guided procedures and interventional applications grow, basic infection control protocols, such as seamless floorings, scrubable surfaces, and hand-washing stations, should be considered.

13. Limits of Applicability and Recommended Design Assistance

The facility design issues identified in this document address only general safety design issues for MRI suites. There are a multitude of site-specific and magnet-specific operational and technical design considerations relevant to MR facility design and construction that are not addressed in these guidelines. These issues include, but are not limited to, patient acuity, staff access, technique conflicts, vibration sensitivity, throughput and efficiency, HIPAA considerations, magnetic contamination, sound transmission, magnet shim tolerances, shielding design, moving metal interferences, MR equipment upgrades, electromagnetic interference, and many others.

In addition to incorporating the guidance from this document, a facility would be well advised to seek expert assistance in the planning and design of MRI and multitechnique radiology suites.

Appendix 3 appears on the next page

APPENDIX 3: Safety Screening Form, MR Hazard Checklist, and Patient Instructions**SAFETY SCREENING FORM FOR MAGNETIC RESONANCE (MR) PROCEDURES**

Date _____
 Name (first middle last) _____
 Female [] Male [] Age _____ Date of Birth _____
 Height _____ Weight _____
 Why are you having this examination (medical problem)? _____

YES NO

_____ Have you ever had an MRI examination before and had a problem?

If yes, please describe: _____

_____ Have you ever had a surgical operation or procedure of any kind?

If yes, list all prior surgeries and approximate dates: _____

_____ Have you ever been injured by a metal object or foreign body (e.g., bullet, BB, shrapnel)?

If yes, please describe: _____

_____ Have you ever had an injury from a metal object in your eye (metal slivers, metal shavings, other metal object)?

If yes, did you seek medical attention? _____

If yes, describe what was found: _____

YES NO

_____ Do you have a history of kidney disease, asthma, or other allergic respiratory disease?

If yes, please describe: _____

_____ Do you have any drug allergies?

If yes, please list drugs: _____

_____ Have you ever received a contrast agent or X-ray dye used for MRI, CT, or other X-ray or study?

If yes, please describe: _____

_____ Have you ever had an X-ray dye or magnetic resonance imaging (MRI) contrast agent allergic reaction?

If yes, please describe: _____

_____ Are you pregnant or suspect you may be pregnant?

_____ Are you breast feeding?

Date of last menstrual period _____


Post-menopausal? _____

Appendix 3 continues on the next page

Safe MR Practices

MR Hazard Checklist

Please mark on the drawing indicating the location of any metal inside your body or site of surgical operation.

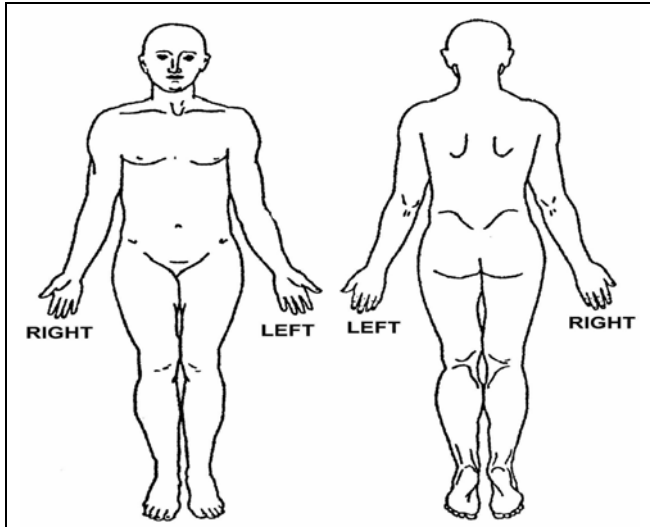

The following items may be harmful to you during your MR scan or may interfere with the MR examination. You must provide a "yes" or "no" for every item. Please indicate if you have or have had any of the following:

YES NO

_____ Any type of electronic, mechanical, or magnetic implant

Type _____

_____ Cardiac pacemaker

_____ Aneurysm clip

_____ Implantable cardiac defibrillator

_____ Neurostimulator

_____ Biostimulator

Type _____

_____ Any type of internal electrodes or wires

_____ Cochlear implant

_____ Hearing aid

_____ Implanted drug pump (e.g., insulin, baclofen, chemotherapy, pain medicine)

_____ Halo vest

_____ Spinal fixation device

_____ Spinal fusion procedure

YES NO

_____ Any type of coil, filter, or stent

Type _____

_____ Any type of metal object (e.g., shrapnel, bullet, BB)

_____ Artificial heart valve

_____ Any type of ear implant

_____ Penile implant

_____ Artificial eye

_____ Eyelid spring

_____ Any type of implant held in place by a magnet

Type _____

_____ Any type of surgical clip or staple

_____ Any IV access port (e.g., Broviac, Port-a-Cath, Hickman, PICC line)

_____ Medication patch (e.g., nitroglycerine, nicotine)

_____ Shunt

_____ Artificial limb or joint

What and where _____

_____ Tissue expander (e.g., breast)

_____ Removable dentures, false teeth, or partial plate

_____ Diaphragm, IUD, pessary

Type _____

_____ Surgical mesh

Location _____

_____ Body piercing

Location _____

_____ Wig, hair implants

_____ Tattoos or tattooed eyeliner

_____ Radiation seeds (e.g., cancer treatment)

_____ Any implanted items (e.g., pins, rods, screws, nails, plates, wires)

_____ Any hair accessories (e.g., bobby pins, barrettes, clips)

_____ Jewelry

_____ Any other type of implanted item

Type _____

Appendix 3 continues on the next page

Instructions for the Patient

1. You are urged to use the ear plugs or headphones that we supply for use during your MRI examination since some patients may find the noise levels unacceptable, and the noise levels may affect your hearing.
2. Remove all jewelry (e.g., necklaces, pins, rings).
3. Remove all hair pins, bobby pins, barrettes, clips, etc.
4. Remove all dentures, false teeth, partial dental plates.
5. Remove hearing aids.
6. Remove eyeglasses.
7. Remove your watch, pager, cell phone, credit and bank cards, and all other cards with a magnetic strip.
8. Remove body piercing objects.

9. Use gown, if provided, or remove all clothing with metal fasteners, zippers, etc.

I attest that the above information is correct to the best of my knowledge. I have read and understand the entire contents of this form, and I have had the opportunity to ask questions regarding the information on this form.

Patient signature _____

MD/RN/RT signature _____

Date _____

Print name of MD, RN, RT _____

For MRI Office Use Only

Patient Name _____

Patient ID Number _____

Referring Physician _____

Procedure _____

Diagnosis _____

Clinical History _____

Hazard Checklist for MRI Personnel**YES NO**

- _____ Endotracheal tube
 _____ Swan-Ganz catheter
 _____ Extraventricular device
 _____ Arterial line transducer

YES NO

- _____ Foley catheter with temperature sensor and/or metal clamp
 _____ Rectal probe
 _____ Esophageal probe
 _____ Tracheotomy tube
 _____ Guidewires

APPENDIX 4: MR Facility Emergency Preparedness Guidelines

Health care facilities have a unique obligation to minimize the disruption from disasters and hasten their ability to restore critical patient care services when interrupted.

Those charged with the operation of MRI facilities have the added complexities of protecting not only the staff and structure, but also the equipment, which may be extraordinarily sensitive to changes in its environment, including vibration, power supply, and water damage.

In the fall of 2005, many watched as Hurricanes Katrina and Rita devastated vast swathes of the U.S. Gulf Coast. Those facilities which were well prepared for the damage, loss of power, and other failures of infrastructure fared far better than those that were not.

Even those not in the likely path of future Gulf hurricanes may have to contend with earthquakes, tornadoes, fires, ice storms, snowstorms, or blackouts, at some point. Particularly those involved in providing patient care should look to how we will provide care at the times when it is most widely and desperately needed. We may find that, while individuals are willing, the facilities, equipment, and infrastructure required to provide clinical care have not been adequately protected.

I. Water Damage

Whether from roof failure, burst pipes, storm surge, or rising rivers, every facility has the potential for water damage to equipment and facilities. Damage can range from inconveniences cured by a couple of hours with a wet-dry vacuum to flooding of equipment electronics. It takes only a small quantity of water in contact with an MRI scanner to incapacitate or destroy the equipment.

To keep leaking roofs, burst pipes, or other overhead damage from dousing MRI equipment, it is recommended that facilities prepare by covering gantries and equipment with sturdy plastic, taped in place, when water damage is an anticipated possibility. To keep processors and gradient cabinets from becoming swamped in a flood situation, electronics that can be lifted off the ground should be moved as far off the floor as possible. RF shields, particularly the floor assembly, may be significantly damaged and need to be replaced in a flood situation if they are not designed to be protected against water damage.

During the 2005 hurricanes, many hospitals and imaging facilities that had emergency generators to help restore power discovered that their sites had generators, or other critical supplies, in basements or other low-lying

areas that were flooded. Facilities should evaluate risks from water damage and assess their preparations for failure of the building enclosure as well as the potential for a flood situation.

2. Structural Damage

MRI presents a particular challenge with structural failure. Although unlikely with current magnet systems, vibrations from seismic events do have the potential to initiate a quench of the magnet system. Structural damage or motion may also damage the RF shield enclosure, potentially degrading image quality until the shield is repaired.

3. Power Outage

Without electrical power to the vacuum pump/cold head to keep the cryogen within a superconducting MRI magnet liquefied, the cryogen will begin to boil off at an accelerated rate. Depending upon cryogen vent design and boil-off rate, the additional cryogenic gas discharge may freeze any accumulated water in the cryogen vent, occluding the pipe and increasing the possibility for a cryogen vent breach in the event of a quench.

At some point, if power to the vacuum pump is not restored, likely a couple days to perhaps a week after power is lost, the magnet will spontaneously quench, discharging most or all of its remaining cryogenic gasses. This poses a safety risk to anyone near the discharge and runs a small but finite risk of potentially permanently damaging the magnet coils.

However, if power to the vacuum pump/cold head and cryogen levels is restored prior to a quench, there should be no long-term consequences to the magnet's operation from a power interruption.

Temporary electrical power may be provided either through on-site or portable generators. Cogeneration, or generating one's own electricity all the time, may not be economically feasible for smaller or stand-alone sites but is increasingly appealing to hospitals for a number of reasons, with emergency capacity being only one.

4. Quench

During the 2005 hurricanes, facilities, fearing extensive damage to their MRI systems from water or protracted power outages, manually initiated preemptive quenches. Under the best circumstances, a quench subjects a magnet to a change of 500°F (260°C) thermal shock within a few dozen seconds, which can cause major physical damage. Rarely, it is possible for the venting cryogenic gases to breach the quench tube and cause significant damage to the magnet room and/or jeopardize the safety of those in the vicinity. At one New Orleans area facility that elected to preemptively quench its magnets, the quench tube reportedly failed and the pressure from the expanding cryogen blew out the control room radiofrequency window (personal communication, Tobias Gilk, October 2005).

Because of the risks to personnel, equipment, and physical facilities, manual magnet quenches are to be initiated only after careful consideration and preparation. In addition to following those specific recommendations provided by the MRI manufacturer, a facility should initiate a preemptive quench in nonemergent situations only after verifying the function of emergency exhaust systems, verifying or providing means of pressure relief, and performing a preliminary visual inspection of the cryogen vent pipe as it leaves the MR unit to check for signs of water or ice inside the pipe (including water leaking from fittings or condensation forming on vent pipe sections).

5. Fire and Police

Though very infrequent, MR suites have been the scene of emergencies requiring fire and/or police response. While it is quite likely this

will be the first time many of the responders have been to an MR suite, this should not be the first time that responding organizations have been introduced to the safety issues for MR. Sites are encouraged to invite police and fire representatives to presentations on MR safety and to provide them with facility tours.

6. Code

In the event that a person within the MR suite should require emergency medical attention, it is imperative that those responding to a call for assistance are aware of, and comply with, MR safety protocols. This includes nurses, physicians, respiratory technicians, paramedics, security personnel, and others.

The impulse to respond immediately must be tempered by an orderly and efficient process to minimize risks to patients, staff, and equipment. This requires specialized training for code teams and, as with fire and police responses, clear lines of authority for screening, access restrictions, and quench authority. Full resuscitation of patients within Zone IV is complicated by the inability to accurately interpret electrocardiographic data. Furthermore, this may place all within Zone IV at risk of injury from ferromagnetic objects which may be on, within, or brought into Zone IV by emergency response personnel responding to a code if one is called in that area. Therefore, after basic cardiopulmonary resuscitation (airway, breathing, chest compressions) is initiated, the patient should be immediately moved out of Zone IV to a prospectively designated location where the code can be run or where the patient will remain until the arrival of emergent response personnel.

It is strongly advised that all MR facilities perform regular drills to rehearse and refine emergency response protocols to protect patients, MR staff, and responders.

7. Prevention

While it is the nature of emergencies to be surprises, we can anticipate the types of incidents that have higher likelihoods given our facilities, practices, and locations. Every facility can anticipate the potential for flooding, fire, and code situations. In addition to these, many areas (e.g., California and coastal Alaska) can expect earthquakes. The central and southern plains states of the United States can anticipate tornados. Colder climates can expect massive snows or ice storms.

State and federal offices of emergency preparedness are dedicated to anticipating and preparing for the specific threats to your region. These offices can serve as an excellent resource regarding risks and strategies for preparation.

Once a disaster has struck, it is important to assess the immediate needs of the community and to restore those critical patient care services first.

Damage to MRI equipment and facilities may not be repaired as quickly. For gravely incapacitated facilities, semitrailer-based MRI units may be the only means of quickly restoring radiology capacity.

All health care facilities should have emergency preparedness plans. The health care plans for MRI facilities should specifically address the unique aspects of MRI equipment. These plans should define who has the authority to authorize nonemergent quenches, procedures for emergency or backup power for the vacuum pump/cold head, as well as instructions on how to protect gantries and sensitive electronics. Facilities should have the necessary supplies pre-positioned and checklists for preparatory and responsive actions. Emergency preparedness plans should also include information necessary for restoring clinical services, including contacts for MRI system vendor, RF

shield vendor, cryogen contractor, MR suite architect and construction contractor, local and state officials, and affiliated hospital and professional organizations.

Below are a few questions that may facilitate the development of an emergency preparedness plan specific to the needs of a facility.

- What are the likely/possible natural disasters to affect the area?
 - What are the likely/possible man-made disasters to affect the area?
 - Is electrical power likely to be interrupted?
 - Would other utilities (natural gas, telecommunications, etc.) likely be interrupted?
 - What equipment would be inoperative during the emergency?
 - What equipment could be damaged by the emergency?
 - What equipment should be provided with critical or backup power?
 - If the utility service is not quickly restored, what other risks may arise?
 - Would patients and staff be able to get to the facility?
 - Would patients or staff be trapped at the facility?
 - How critical is each patient care service provided at the facility?
 - How does the facility protect the equipment needed to support each service?
 - If the facility does not have the resources on site, who can provide them?
-

FOR YOUR INFORMATION

See the accompanying CME/SAM article in this month's issue of *AJR Integrative Imaging* included with your *AJR*.

FOR YOUR INFORMATION

The reader's attention is directed to the commentary on this article, which appears on page 1446.